Efficient, rapid protein purification and on-column cleavage using GSTrap FF columns

Abstract

The need for a purification strategy allowing for the production of soluble, purified proteins in a single step is an invaluable resource and indispensable in many laboratory environments. This Application Note describes a flexible, efficient, and rapid protein purification scheme for the purification of GST-fusion proteins using GSTrap[™] columns and PreScission Protease[™]. The method describes a purification technique that can be applied to reproducibly isolate a variety of diverse proteins of differing function, structure, and chemical nature. This strategy is performed on a single column, yielding "native" proteins (cleaved fusion proteins) in high yield with a high level of enrichment in the 95–97% purity range.

Introduction

The following described purification strategy is an invaluable tool suitable for use in any modern biology laboratory. The method utilizes the technologies of the pGEX-6x vector system, coding for a GST-fusion protein linked with PreScission Protease proteolytic site in conjunction with GSTrap FF affinity columns, and PreScission Protease enzyme. This allows expression and purification of desired target proteins in high yields with high levels of enrichment. Several proteins, which had previously required complex multi-step purification methods, often yielding poorly purified products were selected and assayed. These proteins were sub-cloned, expressed, purified, and cleaved on-column yielding highly pure product. The aim of this strategy is to produce highly pure, "native" protein samples suitable for downstream biochemical assays and structural studies, specifically protein crystallographic analysis.

A flexible method for the purification of diverse gene products

For this method, two unrelated proteins, Pur-x and TLP40, were selected in order to determine the versatility and efficiency of the purification protocol.


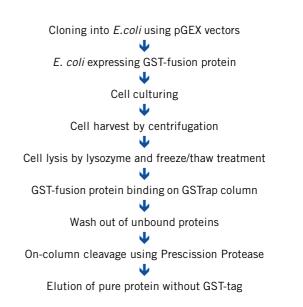


Fig. 1. General purification and on-column cleavage protocol for GST-fusion proteins.

$\mbox{Pur-}\alpha$ is involved in transcription regulation

Pur- α is sequence-specific, single-stranded DNA and RNA binding protein that binds to purine-rich promoter regions with a consensus (GGN)_n sequence (Bergemann and Johnson, 1992). Pur- α has been shown to function in transcriptional regulation by interacting with selective promoter response elements such as: myelin basic protein (MBP) regulatory motif (Haas, S., *et al.*, 1993); neuronal nicotinic acetylcholine receptor (NAChR) genes (Du, Q., *et al.*, 1997); single-stranded cyclic AMP response elementbinding protein (ssCRE-BP; Kuo, C-H., *et al.*, 1999); calmodulin response element (CaMRE) neuropeptide Y (Kuo, C-H., *et al.*, 1999); and human immunodeficiency virus 1 (HIV-1) TAT protein via activation of the TATresponsive DNA element (*upTar*) of the major late promoter

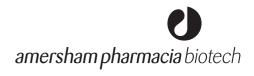

of the JC virus (Krachmarov, C.P., *et al.*, 1996). In addition, Pur- α has been shown to enhance transcription through an RNA element by transactivating HIV-1 through TAR RNA (Chepenik, L.G., *et al.*, 1998).

TLP40 is involved in photosynthesis regulation

TLP40 is an eukaryotic plant protein involved in photosystem II regulation. It is a cyclophillin-like protein with *cis-trans* isomerase activity directed to proline residues, typical for immunophillins (Fulgosi, H., *et al.*, 1998). Biochemical investigations characterizing the physiological role of TLP40 elucidated the interaction with an integral membrane protein phosphatase. This data strongly suggests that TLP40 is a regulator/modulator of phosphorylation states of photosystem II proteins (Vener, A.V., *et al.*, 1999).

Cell growth, protein over-expression, and cell harvesting

The genes encoding for Pur- α and TLP40 protein were independently sub-cloned into the pGEX-6P-1expression vector system and transformed into *E. coli* BL21 over expressing cell lines. Bacteria were grown at +30 °C in LB medium supplemented with 100 µg/ml Carbenicillin in TUNAIR shaking flasks. Cells were grown until OD₆₀₀ = 1.0, and GST-fusion protein synthesis was induced with 1.0 mM IPTG (final concentration), and grown for 3–4 h. Cells were harvested by centrifugation (3000 x g at +4 °C for 30 min), the pellet was resuspended and washed with TB buffer (9.1 mM HEPES, 55 mM MgCl₂, 15 mM CaCl₂, 250 mM KCl adjusted to pH 6.7), centrifuged, and frozen at -80 °C.

Cell lysis

The cell pellet containing GST-fusion protein was resuspended in PBS lysis buffer (140 mM NaCl,

2.7 mM KCl, 10 mM Na₂HPO₄, 1.8 mM KH₂PO₄ adjusted to pH 7.4) supplemented with lysozyme (1 mg/ml, Sigma), CompleteTM Protease Inhibitor Cocktail (1 tablet, Roche), 10 mM MgCl₂, and DNAse I (10 U/ml, Roche). Cells were effectively lysed by repeated (3x) freeze (-170 °C)/thawing (+30 °C). The lysate was centrifuged at 70 000 x g for 25 min at +4 °C. Clarified lysate was then ultracentrifuged at 300 000 x g for 60 min at +4 °C to remove insoluble proteins, cell fragments, and membrane components. The supernatant was applied to a 50 ml SuperloopTM for loading onto a GSTrap FF column.

GST-fusion protein binding

The following described purification procedure was performed on an ÄKTA[™]*explorer* 100 chromatography system. The GST-fusion protein containing supernatant was loaded on two pre-equilibrated GSTrap FF 5 ml columns in series with PBS, pH 7.4 as binding buffer using a Superloop. Sample loading of the column was run at a reduced flow rate (1 ml/min) to optimize the GST affinity media (Glutathione Sepharose[™] 4 Fast Flow) binding interactions on the column. Following GST-fusion protein binding to the column, bound material was washed with PBS, pH 7.4 until the baseline returned. Once the baseline had stabilized, buffer was substituted with PreScission Protease buffer (50 mM Tris HCl, 100 mM NaCl, 1 mM EDTA, 1 mM DTT, pH 8.0).

PreScission buffer was used to more stringently wash the GST-fusion protein bound to the GSTrap FF column, and to pre-equilibrate the column, altering the buffer composition required for efficient PreScission Protease cleavage. PreScission buffer equilibration was continued until both the UV absorbance baseline and the conductivity baseline stabilized, after which the buffer flow was stopped.

On-column GST-fusion protein cleavage with PreScission Protease

With GST-fusion protein effectively bound to the GSTrap FF column, washed and equilibrated in PreScission buffer, the cleavage reaction could be initiated. Exploiting the combined characteristics of the PreScission Protease and the GSTrap FF column, a simple, and rapid on-column GST-fusion protein cleavage can be performed. The reaction proceeds with high efficiency at +4 °C, and is specific and reproducible. The enzyme contains a non-cleavable GST-affinity tag for optimum on-column cleavage, and the column remains on-line, connected to the purification system, eliminating the loss of any material.


Proteolytic digestion requires 2 units enzyme/100 µg of bound GST-fusion protein. PreScission Protease can be diluted in PreScission buffer according to the size of the GSTrap FF column (for example, 4.5 ml for a 5 ml column) and manually injected onto the column at an increased flow rate of 5-7 ml/min. This increased flow rate promotes a reduced affinity interaction time of the PreScission Protease enzyme with the GSTrap FF media allowing for uniformly distributed enzyme throughout the entire column. Following injection, the column is in a closed flow status where efficient PreScission Protease binding occurs, and ideally, even distribution of the enzyme throughout the column is achieved. The system is incubated on-line for 12–16 h at +4 °C.


Elution of cleaved purified protein

Prior to elution, a 1 ml GSTrap FF column (pre-equilibrated with PreScission buffer) was connected downstream to the GSTrap FF proteolytic cleavage column. This tandem column scheme has several functions including acting as a volume delay to route any cleaved material upon flow start-up, minimizing the loss of cleaved product. This allows for a short baseline recalibration before peak elution. The GSTrap FF affinity column also acts as a filter to capture any released cleaved GST protein, uncleaved GST-fusion protein and unbound PreScission Protease.

Cleaved protein elution occurs immediately upon flow startup with PreScission buffer. The elution peak containing the cleaved target protein was eluted first through the 1 ml buffering GSTrap column, then to the detector and fraction collector. The peak profile containing the target protein was somewhat asymmetric in shape due to diffusion effects from the incubation time, and free bound protein. Following target protein elution and return of absorbance baseline, a GST affinity competitor (reduced glutathione) was applied to

> GSTrap 2x5 ml columns in series on ÅKTAexplorer 100 Sample: TLP40 (loaded 50 ml from 2 L culture) Buffer A: PBS, pH 7.4 Buffer B: PreScission buffer, pH 8.0 Buffer C: Reduced glutathione buffer, pH 8.0

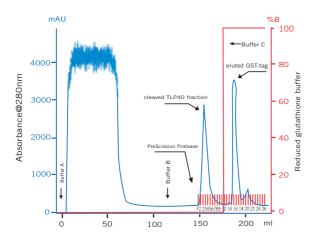

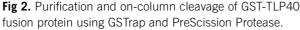

elute GST, unbound GST-fusion protein, and PreScission Protease. A buffer containing reduced glutathione (50 mM Tris-HCl and 10 mM reduced glutathione, pH 8.0) was applied in a one-step gradient (100%) to elute the GST-affinity peak.

The overall purification and on-column cleavage chromatogram can be seen in Figure 2.

Determination of purified target protein

Purification stages and chromatography profiles are analyzed by SDS-PAGE to assay the overall purity and enrichment of the desired target protein, Figure 3. Overall, the final cleaved target protein product is highly enriched with very few contaminants. An estimate of the overall purity is in the range of 95-97% from a single column purification strategy.

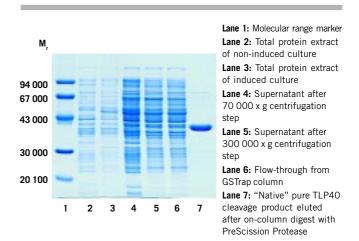


Fig 3. SDS-PAGE analysis of the steps during the purification. Samples were separated on a 3.5–12% polyacrylamide gel. Coomassie[™] staining.

Conclusions

The described protein purification strategy is a simple and rapid method for the isolation of soluble "native" proteins. The GST-fusion protein binding and on-column cleavage assay is greatly facilitated by GSTrap affinity columns and PreScission Protease enzymes allowing for the purification strategy to be performed reproducibly and routinely. The protein yields are highly pure and suitable for downstream processes. This purification strategy exploits the enormous advantages of using the pGEX-6P-1 expression vector system in concert with the engineered PreScission Protease enzyme with its highly specific proteolytic digestion site and uncleavable GST-affinity tag bound to the highly specific binding GSTrap column. With minimal background contaminants, no protease contamination and reproducible results, this method can be regarded as a standard for protein expression and purification trials. A further support for the effectiveness of this strategy invokes a more demanding criterion to judge a protein's purity; that is for it to be able to crystallize. In one of the above test cases, the highly purified protein product crystallizes, reflecting the

0

purification state of the protein, that being highly pure and homogeneous.

This is a method of choice for the purification of "native" proteins due to its simple and rapid strategy with excellent yields.

Acknowledgements: This work was performed by R. Knaust*, S. Eshaghi*, C. Dian*, B. Andersson[#], K. Anzai[†], P. Nordlund* and D. Birse*.

Amersham Pharmacia Biotech thanks specially Professor D. Birse for writing this manuscript and for fruitful discussions.

* The Structural Biochemistry Unit, Department of Biochemistry, Stockholm University, Stockholm, Sweden.

[#] Division of Cell Biology, University of Linköping, Linköping, Sweden.

[†]Department of Biochemistry, College of Pharmacy, Nihon University, Chiba, Japan.

References

- 1. Bergemann, A.D. and Johnson, E.M. (1992) Mol Cell Biol 12; 1257-1265.
- 2. Chepenik, L.G., et al., (1998) Gene 210; 37-44.
- 3. Du, Q. et al., (1997) J Biol Chem 272; 14990-14995.
- 4. Fulgosi, H. et al., (1998) EMBO J 17; 1577-1587.
- 5. Haas, S. et al., (1993) Mol Cell Biol 13; 3103-3112.
- 6. Krachmarov, C.P., et al. (1996) Proc Natl Acad Sci USA 93; 14112-14117.
- 7. Kuo, C.-H., et al. (1999) Biochem Biophys Res Comm 255; 406-411
- 8. Vener, A.V., et al. (1999) Biochemistry 38; 14955-14965.

Ordering information

8		
Designation	Quantity	Code No
GSTrap	2×1 ml	17-5130-02
GSTrap	5×1 ml	17-5130-01
GSTrap	1×5 ml	17-5131-01
Glutathione Sepharose 4		
Fast Flow	25 ml	17-5132-01
Glutathione Sepharose 4		
Fast Flow	100 ml	17-5132-02
Glutathione Sepharose 4		
Fast Flow	500 ml*	17-5132-03
PreScission Protease	500 units	27-0843-01
pGEX-6P-1	25 µg	27-4597-01
pGEX-6P-2	25 µg	27-4598-01
pGEX-6P-3	25 µg	27-4599-01
GST Detection Module	1 kit	27-4590-01
GST 96-Well Detection Module	5 plates	27-4592-01
Anti-GST Antibody	0.5 ml	27-4577-01
Superloop, 50 ml	1	18-1113-82

Larger quantaties available, please contact your Amersham Pharmacia Biotech representative.

(an)

Asia Pacific Tel: +852 2811 8693 Fax: +852 2811 5251 Australasia Tel: +61 2 9894 5152 Fax: +61 2 9899 7511 Austria Tel: 01 576 0616 20 Fax: 01 576 0616 27 Belgium Tel: 0800 73 888 Fax: 03 272 1637 Canada Tel: 1 800 463 5800 Fax: 1 800 567 1008 Central, East, South East Europe Tel: +43 1 982 3826 Fax: +43 1 985 8327 Denmark Tel: 45 16 2400 Fax: 45 16 2424 Finland Tel: 09 512 3940 Fax: 09 512 1710 France Tel: 0169 35 67 00 Fax: 0169 41 9677 Germany Tel: 0761 4903 401 Fax: 0761 4903 405 Haly Tel: 02 27322 1 Fax: 02 27302 212 Japan Tel: 81 3 5331 9336 Fax: 81 5 331 9370 Laxin America Tel: +55 11 3667 5700 Fax: +55 11 3667 87 99 Middle East and Africa Tel: +30 (1) 96 00 687 Fax: +30 (2) 7305 80 410 Fax: 0165 580 401 Fax: 0165 580 401 Norway Tel: 2318 6800 Portugal Tel: 01 417 7035 Fax: 01 417 3184 Russian Federation Tel: +7 (095) 232 0250,956 1137 Fax: +7 (095) 230 6377 South East Asia Tel: 60 3 724 2080 Fax: 018 62 815 494 95 5 Sweden Tel: 018612 19 00 Fax: 018 612 19 10 Switzerland Tel: 01 802 81 50 Fax: 01 802 81 51 UK Tel: 0800 616 928 Fax: 0800 616 927 USA Tel: +1 800 526 3593 Fax: +1 877 295 8102

GSTrap, Sepharose, PreScission, ÄKTA and Superloop are trademarks of Amersham Pharmacia Biotech Limited or its subsidiaries. Amersham is a trademark of Nycomed Amersham plc. Pharmacia and Drop Design are trademarks of Pharmacia & Upjohn. Complete is a trademark of Roche. Coomassie is a trademark of ICl plc. Amersham Pharmacia Biotech AB Björkgatan 30, SE-751 84 Uppsala, Sweden. Amersham Pharmacia Biotech UK Limited Amersham Place, Little Chalfont, Buckinghamshire HP7 9NA, England. Amersham Pharmacia Biotech Inc 800 Centennial Avenue, P0 Box 1327, Piscataway, NJ 08855 USA. Amersham Pharmacia Biotech Limited Amersham Place, Little Chalfont, Buckinghamshire HP7 9NA, England. Amersham Pharmacia Biotech Inc 800 Centennial Avenue, P0 Box 1327, Piscataway, NJ 08855 USA. Amersham Pharmacia Biotech KK. Sanken Building, 3-25-1, Hyakunincho, Shinjuku-ku, Tokyo 169-0073, Japan. All goods and services are sold subject to the terms and conditions of sale of the company within the Amersham Pharmacia Biotech group that supplies them. A copy of these terms and conditions is available on request. © Amersham Pharmacia Biotech A2 000 – All rights reserved.



Produced by Wikströms, Sweden 1001087, 10. 2000 Printed matter. Licence 341 051