



BD TALON™
Metal Affinity Resins
User Manual

PT1320-1 (PR34731)
Published 29 April 2003



Table of Contents

I. Introduction	4
II. List of Components	11
III. Additional Materials Required	13
IV. Buffers for BD TALON™ Purification & Buffer Kits	17
V. Transformation & Protein Expression	18
A. Transformation of Host Cells with Expression Vectors	18
B. Protein Expression	18
VI. Sample Preparation	19
A. BD TALON™ xTractor Buffer Sample Preparation	19
B. Standard Sample Preparation to Isolate Native Proteins	19
C. Standard Sample Preparation to Isolate Denatured Proteins	20
D. Standard Sample Preparation for BD TALON™ CellThru Resin	21
E. Standard HT 96-Well Sample Preparation	22
VII. Protein Purification Protocols	22
A. General Information	22
B. Batch/Gravity-Flow Column Purification BD TALON™ Resin, Superflow Resin, or CellThru Resin	25
C. Large-Scale Batch Purification BD TALON™ Resin, Superflow Resin, or CellThru Resin	26
D. Medium-Pressure Column Purification BD TALON™ Superflow Resin	27
E. BD TALONspin™ Column Purification	28
F. BD TALON™ HT 96-Well Purification Protocol	30
VIII. Resin Washing, Reuse, and Storage	33
IX. Troubleshooting Guide	35
X. References	40
XI. Related Products	41
Appendix A. Reagent Compatibilities and Incompatibilities	42
Appendix B. Mini-Scale Protein Purification Protocol for BD TALON™ or BD TALON™ Superflow Resin	43
Appendix C. Vector Information	45

Table of Contents *continued*

List of Figures

Figure 1. Schematic diagram of the BD TALON™ IMAC System.	5
Figure 2. Elution mechanism of recombinant polyhistidine-tagged proteins from BD TALON™ Resin.	6
Figure 3. Binding of histidines to the BD TALON™ Resin metal ion.	6
Figure 4. Purification of polyhistidine-tagged proteins	15
Figure 5. pHAT10/11/12 Combined Vector Map and MCS.	45
Figure 6. pHAT20 combined vector map and MCS.	46

List of Tables

Table I. Protein purification using BD TALON™ Resins.	8
Table II. BD TALON Resin characteristics	10
Table III. Reagent Compatibility	42

Notice to Purchaser

This product is intended to be used for research purposes only. It is not to be used for drug or diagnostic purposes nor is it intended for human use. BD Biosciences Clontech products may not be resold, modified for resale, or used to manufacture commercial products without written approval of BD Biosciences Clontech.

BD TALON™ products are covered under U.S. Patent No. 5,962,641.

Sepharose® is a registered trademark of Pharmacia LKB Biotechnology.

Triton® is a registered trademark of Rohm and Haas Co.

Superflow™, Uniflow™ and CellThru™ are trademarks of Sterogene Bioseparations, Inc.

BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. ©2003 BD

I. Introduction

Proteins have evolved very complex structures in order to perform a diverse array of functions. As a result, their physicochemical properties vary greatly, posing difficulties for developing versatile purification protocols. One way to circumvent this problem is to incorporate a purification tag into the primary amino acid sequence of a target protein, thus constructing a recombinant protein with a binding site that allows purification under well-defined, generic conditions.

Immobilized Metal Affinity Chromatography (IMAC)

IMAC was introduced in 1975 as a group-specific affinity technique for separating proteins (Porath *et al.*, 1975). The principle is based on the reversible interaction between various amino acid side chains and immobilized metal ions. Depending on the immobilized metal ion, different side chains can be involved in the adsorption process. Most notably, histidine, cysteine, and tryptophan side chains have been implicated in protein binding to immobilized transition metal ions and zinc (Figure 1, Porath, 1985; Sulkowski, 1985; Hemdan & Porath, 1985a; Hemdan & Porath, 1985b; Zhao *et al.*, 1991).

BD TALON™ IMAC Resins

BD TALON™ Resins are durable, cobalt-based IMAC resins designed to purify recombinant polyhistidine-tagged proteins (Bush *et al.*, 1991). These resins are compatible with many commonly used reagents (Appendix A), and allow protein purification under native or denaturing conditions. They can be used with all prokaryotic and eukaryotic expression systems in a variety of formats, including small- (or mini-) scale batch screening, large-scale batch preparations, and methods using gravity-flow columns and spin columns. In addition, protocols used with Ni⁺²-based IMAC columns usually work with BD TALON™ resins.

Tetradentate metal chelator

To overcome the problem of metal leakage encountered with other IMAC resins, BD TALON™ Resin utilizes a special tetradentate metal chelator for purifying recombinant polyhistidine-tagged proteins (U.S. Patent No. 5,962,641). This chelator tightly holds the electropositive metal in an electronegative pocket (Figure 1), which is ideal for binding metal ions such as cobalt. The binding pocket is an octahedral structure in which four of the six metal coordination sites are occupied by the BD TALON Resin ligand. This process enhances the protein binding capacity of BD TALON Resin by making the bound metal ion accessible to surrounding polyhistidine-tagged proteins. The tetradentate metal binding means that no metal loss occurs during protein purification under recommended conditions, even in the presence of strong denaturants such as 6 M guanidinium. Such durability allows BD TALON Resin to be reused (See Section VIII).

Cobalt IMAC Resin permits milder elution conditions

BD TALON Resin exhibits subtle yet important differences in binding of polyhistidine-tagged proteins when compared with nickel IMAC resins. For example, nickel-based IMAC resins often exhibit an undesirable tendency to bind

I. Introduction *continued*

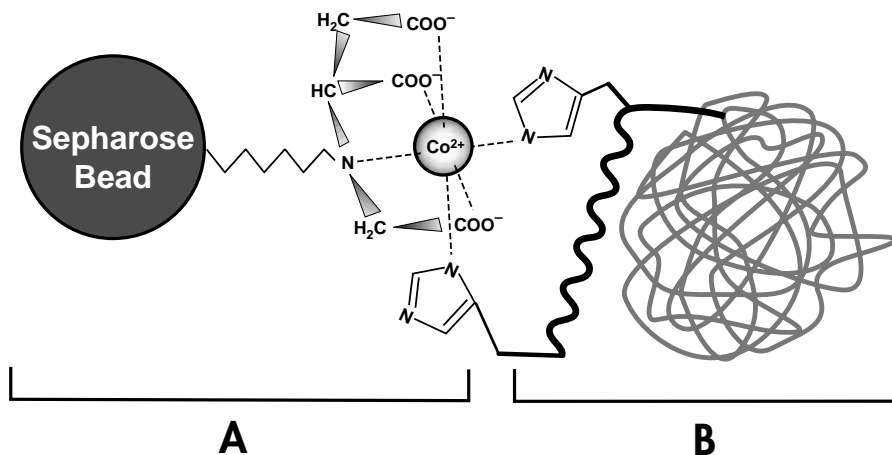


Figure 1. Schematic diagram of the BD TALON™ IMAC System. Part A. BD TALON Metal Affinity Resin; A Sephadex bead bearing the tetradentate chelator of the Co²⁺ metal ion. **Part B.** The polyhistidine-tagged recombinant protein binds to the resin.

unwanted host proteins containing exposed histidine residues (Kasher *et al.*, 1993). While BD TALON Resin binds polyhistidine-tagged proteins with enhanced selectivity over nickel-based resins, it exhibits a significantly reduced affinity for host proteins. This behavior offers two practical advantages. First, virtually no background proteins are bound to the resin when the sample is applied; consequently, cumbersome washing procedures are not generally required before protein elution. Second, polyhistidine-tagged proteins elute from the resin under slightly less stringent conditions—a slightly higher pH or lower imidazole concentration—than with nickel IMAC resins. Elution occurs when the imidazole nitrogen (pKa of 5.97) is protonated (Figure 2), generating a positively charged ammonium ion, which is repelled by the positively charged metal atom. Alternatively, the bound polyhistidine-tagged protein can be competitively eluted by simply adding imidazole to the elution buffer, because imidazole is identical to the histidine side chain.

Polyhistidine affinity tags

Histidines exhibit highly selective coordination with certain transition metals and have great utility in IMAC. Under conditions of physiological pH, histidine binds by sharing electron density of the imidazole nitrogen with the electron-deficient orbitals of transition metals (Figure 3). Although three histidines may bind transition metals under certain conditions, six histidines reliably bind transition metals in the presence of strong denaturants such as guanidinium (Hochuli *et al.*, 1987). Such protein tags are commonly referred to as “6xHis,” “hexaHis,” or “(His)₆.”

I. Introduction *continued*

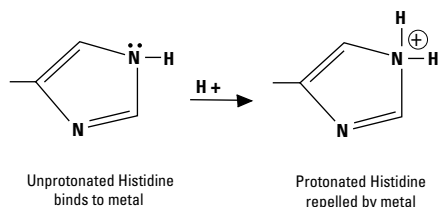


Figure 2. Elution mechanism of recombinant polyhistidine-tagged proteins from BD TALON™ Resin. Elution occurs when the imidazole nitrogen (pKa = 5.97) is protonated, generating a positively charged ammonium ion which is repelled by the positively charged metal ion. Alternatively, the bound polyhistidine-tagged protein can be competitively eluted by adding imidazole to the elution buffer.

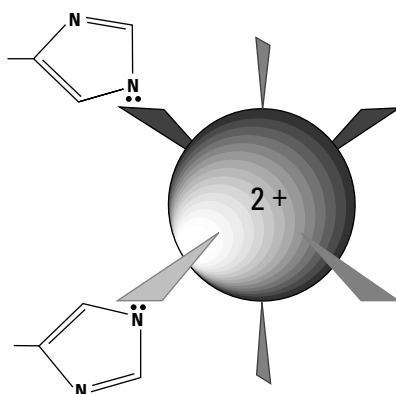


Figure 3. Binding of histidines to the BD TALON™ Resin metal ion. Under conditions of physiological pH, histidine binds by sharing imidazole nitrogen electron density with the electron-deficient orbitals of the metal ion.

HAT—a novel IMAC affinity tag

With the advent of recombinant genetic technologies, the design and production of recombinant proteins containing novel polyhistidine tags on their N- or C-termini has become more straightforward (Hochuli *et al.*, 1987; Hochuli *et al.*, 1988). The HAT sequence (patent pending) is a novel IMAC affinity tag derived from a unique natural protein sequence (Chaga *et al.*, 1999). It contains six histidines unevenly interleaved by other amino acid residues (see Appendix C). The HAT amino acid sequence is derived from the N-terminus of chicken muscle lactate dehydrogenase—a sequence that is unique among reported protein sequences. The novel tag does not have the excessive positive charge characteristic of the commonly used 6xHis tag, thus contributing to better solubility of HAT-fusion proteins and similar affinity towards immobilized transition metal ions and zinc. BD Biosciences Clontech offers the HAT Protein Expression and Purification System (Cat. No. 631205)—a complete system containing reagents and vectors designed for bacterial expression and purification of HAT (histidine affinity tag) proteins. Each of the three vectors—pHAT10, pHAT11, and pHAT12—contain a multiple cloning site (MCS) in all three frames to allow cloning of target cDNA. (For vector map and MCS, see Appendix C of this User Manual.) A conveniently located enterokinase proteolytic site between the HAT sequence and the MCS provides a means for removing the affinity tag.

For more information, see the HAT Expression & Purification System User Manual (PT3250-1), which can be downloaded from our web site at www.bdbiosciences.com/clontech/.

I. Introduction *continued*

Overview of BD TALON™ Resins

The following is a list of different resin formats to meet your purification needs.

- **BD TALON™ Metal Affinity Resin** is useful for batch and low-pressure chromatographic applications. This resin utilizes Sepharose® CL-6B (Pharmacia LKB Biotechnology), a durable substrate that performs very well under native and denaturing conditions in centrifuge-mediated purification schemes. The large pore size resin has a high-binding capacity. This resin is also available pre-packed in 2-ml gravity columns or in a high-throughput (HT) 96-well plate (See List of Components for more information).

- **BD TALON™ Superflow Resin** is useful for a range of applications, including medium pressure applications with FPLC systems at back pressures of up to 150 psi (1 MPa) and high flow rates up to 5 ml per cm² per min. This resin is recommended if short purification times are essential, or if purification protocols developed at bench scale will be scaled up for larger volumes.

This resin utilizes Superflow-6 (Sterogene Bioseparations, Inc.), an agarose-based medium featuring a unique polysaccharide composition that resists biological degradation. Superflow-6 beads are also stabilized by a chemical cross-linking reaction that allows flow rates up to 10 times higher than are possible with regular cross-linked beads.

- **BD TALON™ CellThru Resin** is a novel IMAC resin for purifying polyhistidine-tagged proteins from crude cell lysates, sonicates, and fermentation liquids. The larger bead size of BD TALON™ CellThru Resin (300–500 μm) permits cellular debris to flow through the column, eliminating the need for high-speed centrifugation. Destabilizing factors are removed more quickly with this resin than with other IMAC resins, because the number of steps are reduced.

CellThru 2-ml & 10-ml Disposable Columns have a large filter pore size (90–130 μm) that allows cellular debris to flow through the column during the purification process. The 2-ml columns are suitable for 1–2 ml bed volumes, while the 10-ml columns are suitable for 5–10 ml bed volumes.

- **BD TALONspin™ Columns** are ideal for rapidly and simultaneously purifying small amounts of polyhistidine-tagged proteins. These columns are recommended for single-use applications or for use as mini gravity-flow columns. Each column contains 0.5 ml of BD TALON-NX™ Resin, which is optimized for performance in a spin column. Each column will yield 2–4 mg of polyhistidine-tagged protein; exact yields will vary with conditions used and polyhistidine-tagged protein characteristics. In addition, yield and purity will depend upon expression level and lysate concentration. Beginning with the clarified sample, the entire procedure takes approximately 30 minutes.

I. Introduction *continued*

TABLE I. PROTEIN PURIFICATION USING BD TALON RESINS

Method	Application	Key Benefit
BD TALON™ Metal Affinity Resin or BD TALON™ Superflow Resin		
Mini-Scale (Appendix B)	<ul style="list-style-type: none"> • Check for presence of tagged protein • Estimate expression levels • Test buffer conditions 	<ul style="list-style-type: none"> • Fast • Requires only 1 ml of cell culture + 1 ml of resin
Batch/Gravity Flow Column (Sec. VII.B)	<ul style="list-style-type: none"> • Purify ≥ 5 mg of tagged protein using 1 ml of resin 	<ul style="list-style-type: none"> • Very high purity • Does not require pressurized column equipment
Large-Scale (Sec. VII.C & D)	<ul style="list-style-type: none"> • Large- and production-scale purification; easy to scale up 	<ul style="list-style-type: none"> • Faster than protocols that use gravity-flow columns • Higher purity than using batch process alone
BD TALON™ CellThru Resin		
Batch/Gravity Flow Column & Large-Scale (Sec. VII.B & C)	<ul style="list-style-type: none"> • For purifying proteins from nonclarified cell lysates, sonicates, or fermentation liquids 	<ul style="list-style-type: none"> • Fast • Does not require high-speed centrifugation
BD TALONspin™ Columns		
Spin Column (Sec. VII.E)	<ul style="list-style-type: none"> • Process several different samples simultaneously • Obtain 2–4 mg of purified protein per spin column 	<ul style="list-style-type: none"> • Fast (~30 min)* • Uses only 0.6–1 ml of cell culture
BD TALON™ HT 96-Well Plates		
96-well Plates (Sec. VII.F)	<ul style="list-style-type: none"> • High-throughput processing of samples • Obtain up to 1.0 mg of purified protein per well 	<ul style="list-style-type: none"> • Fast (<30 min)* • Uses up to 2 ml of crude lysate per load

* Starting with clarified lysate; does not include time to analyze samples.

I. Introduction *continued*

Protein Purification Methods Using BD TALON™ Resin

The following general guidelines are used for purifying polyhistidine-tagged protein from transformed *E. coli* cultures. Table I provides an overview of BD TALON™ Resin protein purification methods and applications. Choose a method that best suits your research needs.

- Use 2 ml of resin suspension per ~3 mg of anticipated polyhistidine-tagged protein. 2 ml of homogeneously resuspended resin will provide 1 ml (bed volume) of BD TALON Resin.
- The buffers and purification conditions should work well for most soluble, monomeric proteins expressed in *E. coli*.
- Initially, test each different expression system and polyhistidine-tagged protein in small-scale batch purification to determine expression levels and to optimize the protocol. A mini-scale batch purification protocol is provided in Appendix B; alternatively, you can use a BD TALONspin Column.
- Purification methods that work with Nickel or Zinc-based IMAC resins should also work with these resins. However, some optimization may be required.

Note: BD TALON resin has been optimized and should only be used with the buffer formulations outlined in this user manual for optimal performance.

Choosing the Buffers: Imidazole Versus pH Gradient

BD TALON™ Resin purification schemes typically use either an imidazole or a pH gradient for washing and elution. Imidazole in the Equilibration and/or Equilibration/Wash Buffers minimizes nonspecific binding and reduces the amount of contaminating proteins. Thus, we recommend first purifying polyhistidine-tagged proteins using an imidazole gradient. However, imidazole and polyhistidine-tagged proteins absorb at 280 nm and elution peaks may be difficult to detect spectrophotometrically, especially if you are purifying small amounts of polyhistidine-tagged proteins. In these cases, collect the leading edge of the imidazole breakthrough peak and check for polyhistidine-tagged proteins by a protein specific assay (Bradford, 1976) and SDS/PAGE. Alternatively, use a pH gradient to purify polyhistidine-tagged proteins that are stable from pH range 5.0–7.0. See Section III for buffer compositions.

Elution Strategy: Step Versus Linear Gradients

In most cases, step gradients are preferred over linear gradients, because linear gradients lead to broad elution peaks, which can dilute the product and make detection difficult. Scaling-up step gradients is also less complicated than scaling-up linear gradients.

I. Introduction *continued*

TABLE II. BD TALON™ RESIN CHARACTERISTICS

Feature	BD TALON™ Resin	BD TALON™ Superflow	BD TALON™ CellThru	BD TALONspin™ Column
Capacity (mg protein/ml resin)	5–10	5–14	5–10	2–4
Matrix	Sepharose CL-6B	Superflow	Uniflow	Silica
Bead size (µm)	45–165	60–160	300–500	16–24
Max. Linear flow rate (cm/hr)	75–150	3,000	800	n/a
Max. Vol. flow rate* (ml/min)	0.5–1.0	50	13	n/a
Max. Pressure	2.8 psi 0.2 bar 0.02 MPa	140 psi 10 bar 0.97 MPa	9 psi 0.62 bar 0.02 MPa	n/a
pH stability	2–14 (2 hr) 3–14 (24 hr)	2–14 (2 hr) 3–14 (24 hr)	2–14 (2 hr) 3–14 (24 hr)	2–8.5 (2 hr) 2–7.5 (24 hr)
Protein exclusion limit (Da)	4 x 10 ⁷	4 x 10 ⁶	2 x 10 ⁷	n/a

*determined on a 5 x 1 cm column

n/a = not applicable

II. List of Components

BD TALON Resin, BD TALON Superflow Resin, and BD TALON CellThru Resin are supplied as 50% (w/v) slurries in nonbuffered 20% ethanol. Please note that during shipping and storage, the resin will settle; thus, we recommend that you thoroughly resuspend it before aliquotting. 2 ml of homogeneously resuspended resin will provide 1 ml of BD TALON Resin with a binding capacity of at least 5 mg of polyhistidine-tagged protein.

Store all of these resins, columns and buffers at 4°C unless otherwise indicated.
Do not freeze BD TALON™ Resins.

- **BD TALON™ Metal Affinity Resin**

<u>Cat. No.</u>	<u>Amount</u>
-----------------	---------------

635501	10 ml
635502	25 ml
635503	100 ml
635504	250 ml

- **BD TALON™ Superflow Resin**

<u>Cat. No.</u>	<u>Amount</u>
-----------------	---------------

635506	25 ml
635507	100 ml

- **BD TALON™ HT 96-Well Plate** (Cat. No. 635622)

1	BD TALON 96-Well Plate
1	Plate Top Seal
1	Plate Base Seal
1	Collection Plate

- **BD TALONspin™ Columns** (Cat. Nos. 635601, 635602, 635603, 635604)

These columns contain 0.5 ml of BD TALON-NX resin as a 50% suspension in nonbuffered 20% ethanol.

- **BD TALON™ CellThru Resin**

<u>Cat. No.</u>	<u>Amount</u>
-----------------	---------------

635509	10 ml
635510	100 ml

- **BD TALON™ CellThru Disposable Columns**

<u>Cat. No.</u>	<u>Size</u>
-----------------	-------------

635512	2-ml column
635513	10-ml column

- **BD TALON™ 2-ml Disposable Gravity Columns** (Cat. No. 635606)

II. List of Components *continued*

- **BD TALON™ Purification Kit** (Cat. No. 635515)

- 10 ml BD TALON™ Metal Affinity Resin
- 160 ml 5X Equilibration/Wash Buffer
(250 mM Sodium Phosphate, 1.5 M Sodium Chloride, pH 7)
- 160 ml 5X Equilibration Buffer
(250 mM Sodium Phosphate, 1.5 M Sodium Chloride, pH 8)
- 25 ml 10X Elution Buffer
(1.5 M Imidazole, pH 7)
- 5 2-ml Disposable Gravity Columns
- 1 10-ml Disposable Gravity Column

- **BD TALON™ Buffer Kit** (Cat. No. 635514)

- 160 ml 5X Equilibration/Wash Buffer
(250 mM Sodium Phosphate, 1.5 M Sodium Chloride, pH 7)
- 160 ml 5X Equilibration Buffer
(250 mM Sodium Phosphate, 1.5 M Sodium Chloride, pH 8)
- 25 ml 10X Elution Buffer
(1.5 M Imidazole, pH 7)

- **BD TALON™ xTractor Buffer Kit** (Cat. No. 635623)

Store DNase I at -20°C .

If a precipitate has formed in the Lysozyme solution, allow the tube to warm at room temperature and gently invert the tube. The solution may remain turbid after this procedure.

- 200 ml xTractor Buffer
- 2.5 ml 50X Lysozyme
- 400 μl DNase (1 U/ μl)

- **BD TALON™ xTractor Buffer** (Cat. No. 635625)

- 500 ml xTractor Buffer

III. Additional Materials Required

See Section IV for preparing buffers with the BD TALON™ Purification Kit (Cat. No. 635515) or the BD TALON™ Buffer Kit (Cat. No. 635514). If you have **not purchased** those kits, we recommend preparing the following buffers for purifying polyhistidine-tagged proteins under native or denaturing conditions. Before preparing other buffer compositions, please consult Appendix A to evaluate resin compatibility.

Choosing Buffers

To decrease the amount of nonspecifically bound proteins, we recommend using the **Equilibration/Wash Buffer** at pH 7.0 during purification; however, if your target protein is more stable at pH 8.0, or if it does not adsorb at pH 7.0, use the Equilibration Buffer at pH 8.0 (in place of the Equilibration/Wash Buffer) during all extraction and wash steps. Note that at elevated pH values, amino acids other than histidine, as well as the peptide bond, contribute to protein adsorption. Thus, proteins without a polyhistidine tag can also adsorb to IMAC Resins, which decreases resin capacity and the final purity of your target protein. You may choose to use either native or denaturing buffer conditions, depending on the solubility of your protein. Figure 4 outlines the purification procedure.

A. Native Buffers

Native protein purification regimens use buffer conditions that preserve the native, three-dimensional structure and surface charge characteristics of a selected soluble protein during harvest from an expression host. The low affinity of BD TALON™ Resin for nonpolyhistidine-tagged proteins minimizes contaminant carryover. In addition, increasing buffer ionic strength can minimize nonspecific interactions. Regardless of the conditions used and the nature of the polyhistidine-tagged protein being purified, most applications will benefit from the presence of 100–500 mM NaCl in the IMAC buffer. In many cases, adding glycerol or ethylene glycol neutralizes nonspecific hydrophobic interactions. Small amounts of nonionic detergent may also dissociate weakly bound species.

- **Phosphate Buffer (PB; pH 7.5)**

	<u>Final conc.</u>	<u>To prepare 2 L of solution</u>
Na ₂ HPO ₄	58 mM	16.5 g
NaH ₂ PO ₄	17 mM	4.1 g
NaCl	68 mM	8.0 g

Dissolve the above components in 1.8 L of deionized H₂O. Adjust to pH 7.5 with 0.1 N NaOH. Add deionized H₂O to final volume of 2 L. Store at room temperature.

III. Additional Materials Required *continued*

- **1X Equilibration/Wash Buffer** (pH 7.0)
 - 50 mM Sodium Phosphate
 - 300 mM NaCl
- **1X Equilibration Buffer** (pH 8.0)
 - 50 mM Sodium Phosphate
 - 300 mM NaCl
- **1X Elution Buffer** (See also Section VI.B)
 - **Imidazole Elution** (pH 7.0)
 - 50 mM Sodium Phosphate
 - 300 mM NaCl
 - 150 mM Imidazole
 - **pH Elution** (pH 5.0)
 - 50 mM Sodium Acetate
 - 300 mM NaCl
- **HT 96-Well Plate Wash Buffer** (pH 7.0)
 - 83 mM Sodium Phosphate
 - 500 mM NaCl
 - 10 mM Imidazole

B. Denaturing Buffers

Denaturants, such as 6 M guanidinium, enhance protein solubility. Because proteins overexpressed in prokaryotic systems are sometimes insoluble, you may need to purify proteins under denaturing conditions. When purifying proteins under denaturing conditions, we recommend preparing the buffers indicated below.

In the presence of 6 M guanidinium, *the resin's color will change from a pinkish-mauve to violet due to a change in metal ion hydration in response to the chaotrope*. After removal of the chaotrope, the resin will return to a pinkish-mauve color. The change to violet does not reflect any change in the physical or chemical properties of the resin. In fact, the color change can be useful for indicating the buffer in which the resin is suspended, and for following the movement of guanidinium through the resin bed.

- **1X Equilibration/Wash Buffer** (pH 7.0)
 - 50 mM Sodium Phosphate
 - 6 M Guanidine-HCl
 - 300 mM NaCl
- **1X Equilibration Buffer** (pH 8.0)
 - 50 mM Sodium Phosphate
 - 6 M Guanidine-HCl
 - 300 mM NaCl

III. Additional Materials Required *continued*

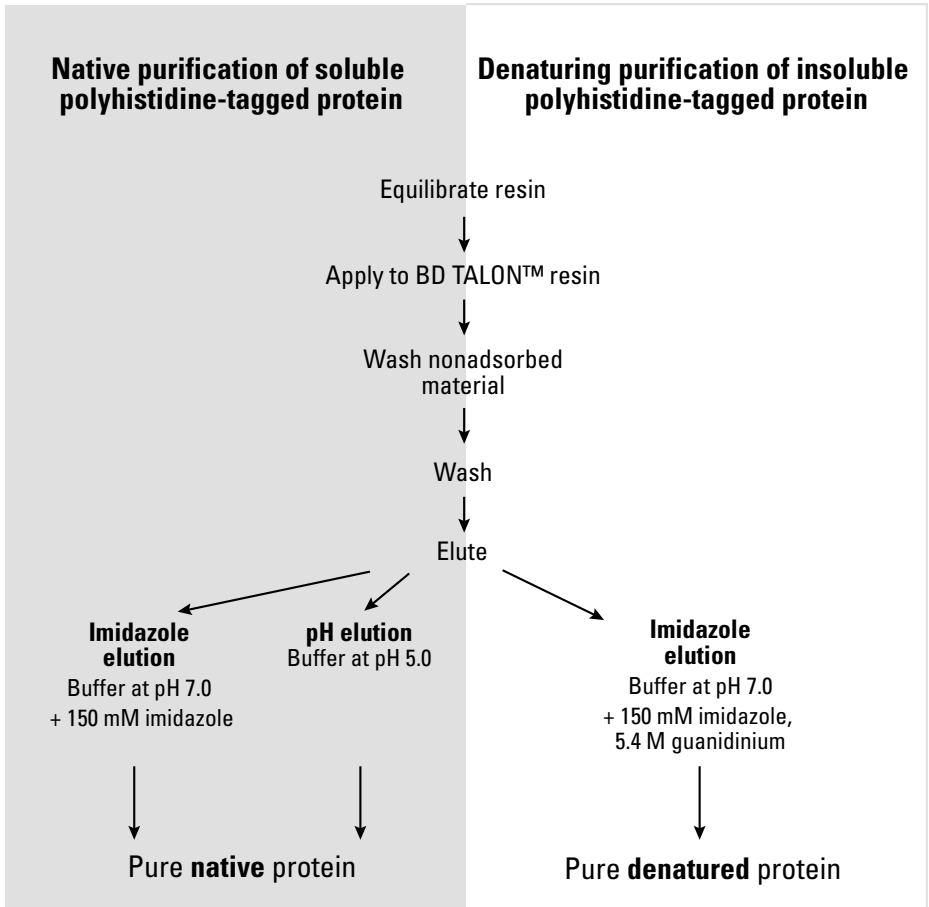


Figure 4. Purification of polyhistidine-tagged proteins using BD TALON™ Resin. The protocols in this User Manual are designed using the Equilibration/Wash Buffer at pH 7.0. If your target protein is more stable at pH 8.0, or if it does not adsorb at pH 7.0, use the Equilibration Buffer (pH 8.0) instead of the Equilibration/Wash Buffer during the extraction and wash steps.

- **1X Imidazole Elution Buffer (pH 7.0)**
 - 45 mM Sodium Phosphate
 - 5.4 M Guanidine-HCl
 - 270 mM NaCl
 - 150 mM Imidazole

III. Additional Materials Required *continued*

C. Additional Buffers & Reagents

- **MES Buffer**

20 mM 2-(N-morpholine)-ethanesulfonic acid (MES), pH 5.0

- **5X SDS PAGE Sample Buffer**

15% β -Mercaptoethanol (β -ME)

15% SDS

50% Glycerol

1.5% Bromophenol blue

- **Imidazole** (Sigma, Cat. No. I0250) Also suitable for FPLC Applications

- **Bio-Rad Protein Assay** (Bio-Rad, Cat. No. 500-0001)

D. Additional Materials required for BD TALON™ CellThru

- CellThru 2-ml Disposable Columns (Cat. No. 635512)

- CellThru 10-ml Disposable Columns (Cat. No. 635513)

E. Additional Materials for BD TALON™ HT 96 Plate

Vacuum purification

- NucleoVac (Cat. No. 636030), QIAVac 96 (Qiagen Cat. No. 19504) or similar vacuum manifold

- (Extra) Collection 96-Deep Well Titer Plates (Whatman Cat. No. 7701-5200 or Evergreen Cat. No. 240-8556-030)

Centrifugation

- Centrifuge with a rotor for centrifugation of microtiter plates, such as the Allegra 6R Centrifuge (Beckman Coulter) with the GH 3.8; GH 3.8A; or JS 4 Beckman rotors.

IV. Buffers for BD TALON™ Purification and Buffer Kits

If you have purchased the **BD TALON™ Purification or Buffer Kits**, prepare buffers as described below. To decrease the amount of nonspecifically bound proteins, we recommend using the **Equilibration/Wash Buffer** at pH 7.0 during purification; however, if your target protein is more stable at pH 8.0, or if it does not adsorb to the resin at pH 7.0, use the Equilibration Buffer (pH 8.0) in place of the Equilibration/Wash Buffer during all extraction and wash steps. Note that at elevated pH values, amino acids other than histidine, as well as the peptide bond, can be adsorbed by BD TALON Resins; Thus, in high pH conditions (pH>8.0) proteins without a polyhistidine tag can be adsorbed and decrease resin capacity and the final purity of your target protein.

Note: If precipitate is observed in the buffers, warm them at 37°C, and stir or shake to dissolve precipitate prior to diluting and using the buffers.

A. BD TALON xTractor Buffer: No preparation necessary except optional addition of DNase I or Lysozyme (see Section VI.A).

B. Equilibration Buffers

1. Dilute one part of the 5X Equilibration/Wash Buffer or 5X Equilibration Buffer with four parts of deionized water.
2. Check and correct pH if necessary. The 1X Equilibration/Wash Buffer should be pH 7.0, while the 1X Equilibration Buffer should be pH 8.0.

C. Elution Buffer

Dilute one part of the 10X Elution Buffer with nine parts of 1X Equilibration/Wash Buffer (pH 7.0) (or 1X Equilibration Buffer [pH 8.0], depending on the solubility of your protein) prepared in Step A.

D. Denaturing Conditions

Add 6 M Guanidinium to the Equilibration/Wash Buffer (pH 7.0), or Equilibration Buffer (pH 8.0), and the Elution Buffer prepared in Steps A and B, respectively.

Note: Perform all steps during the purification procedure in the presence of 6 M Guanidinium. Protein samples containing high guanidinium concentrations form a precipitate when loaded on SDS polyacrylamide gels. Therefore, dialyze the sample overnight in a buffered solution containing 8 M urea before loading it onto the gel.

E. Wash Buffers

- In general, use the Equilibration/Wash Buffer at pH 7.0 to wash nonadsorbed proteins. If the protein is not stable at pH 7.0, then use the Equilibration Buffer at pH 8.0 with 5–10 mM imidazole.
- If your host cell line produces unwanted multi-histidine proteins, incorporate a more stringent wash:
Dilute 5X Elution Buffer in either 1X Equilibration/Wash Buffer or 1X Equilibration Buffer for a final concentration of 5–10 mM Imidazole (1:300–1:150).
- For HT 96-Well Plate Wash Buffer: Dilute Elution Buffer 1:150 in 3X Equilibration/Wash Buffer.

V. Transformation & Protein Expression

A. Transformation of Host Cells with Expression Vectors

The following protocol is for chemically-induced transformation of *E. coli* competent cells. Perform control transformations in parallel.

Note: Use JM109 or another lac-inducible cell line to see induction of expression. For tighter control of expression levels, use our PROTeT 6xHN Bacterial Expression System—especially recommended for expression of cytotoxic proteins.

1. On ice, thaw a tube containing 100 μ l of 0.5 M β -mercaptoethanol (β -ME) and one 50- μ l tube of frozen *E. coli* competent cells for each ligation/transformation.
2. Dispense 2 μ l of 0.5 M β -ME into each tube of competent cells and mix.
3. Dispense 2 μ l of DNA directly into the mixture from Step 2.
4. Incubate tubes on ice for 30 min.
5. Heat shock for exactly 30 sec in a 42°C water bath.
6. Remove tubes from water bath and place on ice for 2 min.
7. Add 250 μ l of SOC medium to each tube at room temperature.
8. Shake the tubes horizontally at 37°C for 1 hr at 225 rpm in a rotary shaking incubator.
9. Spread transformation mixtures onto LB-ampicillin (50 μ g/ml) agar plates [containing X-gal (75 μ g/ml) and IPTG (1 mM) if blue-white selection is desired]. Incubate the plates at 37°C overnight.

B. Protein Expression

1. Grow an overnight culture of *E. coli* transformed with your expression plasmid. If you can isolate a sufficient amount of protein from this culture, proceed to Step 3 after taking a 1-ml sample for electrophoretic analysis. Centrifuge the sample at 1,000–3,000 \times g for 15 min at 4°C, remove the supernatant, and store the cell pellet at –20°C.

Note: If a large-scale preparation of the protein is required, proceed to Step 2.

2. If you need a greater quantity of the target protein, use 20 ml of overnight culture to inoculate 1 L of medium. Incubate with shaking for another 1–2 hr, until the culture has an absorbance of \sim 0.6 OD₆₀₀. Remove a 1-ml sample of the culture, centrifuge at 1,000–3,000 \times g for 15 min at 4°C, remove the supernatant and store the cell pellet at –20°C for electrophoretic analysis.
3. Induce expression by adding an appropriate inducer. For example, the *lac* promoter in the pHAT10 expression vector can be induced with 1 mM IPTG. Continue the incubation for another 3–5 hr.
4. Remove a 1-ml sample of the culture, centrifuge at 1,000–3,000 \times g for 15 min at 4°C, remove the supernatant, and store the cell pellet at –20°C for electrophoretic analysis.
5. Proceed with sample preparation (Section VI).

VI. Sample Preparation

A. BD TALON™ xTractor Buffer Sample Preparation

This procedure has been optimized for extraction of native proteins from fresh or frozen bacterial cell pellet using BD TALON xTractor Buffer. The volumes of this extraction can be adjusted, as long as 20 ml of xTractor Buffer are used per 1 g of cell pellet.

1. Add 20 ml of xTractor Buffer to 1 g of bacterial cell pellet. Mix gently. Pipet the mixture up and down to fully resuspend the pellet.

Notes

- For BD TALON HT 96-Well Plate (Cat. No. K1254-1), resuspend 40–100 mg of bacterial cell pellet in 2 ml of xTractor Buffer.
 - For a ~1 g bacterial pellet, we harvest approximately 500 ml of *E. coli* culture, which had been induced at log phase (O.D.=0.6–0.8) and then incubated an additional 2–4 hr for protein production.
2. [Optional]: Add 40 μ l of 1U/ μ l DNase I solution and 200 μ l of 50X Lysozyme solution.

Notes

- DNase I reduces the viscosity of the lysate allowing for more efficient removal of cellular debris. DNase I can be used without lysozyme. However, if you are treating cells with lysozyme, then you **must** treat cells with DNase I as well.
 - Lysozyme helps to fully disrupt bacterial walls, and thus it has been demonstrated to be highly beneficial in extraction of high molecular weight proteins (>40 kDa). However, lysozyme should be omitted for mammalian extraction procedures as well as when lysozyme interferes with your protein's functionality.
3. Mix gently, pipetting up and down several times.
 4. Incubate with gentle shaking or stirring for 10 min at room temperature. (You may incubate the solution at 4°C, if desired).

Note: At the end of this incubation period, there should be no visible particles. If cell pellet fragments are present, resuspend them by pipetting solution up and down and incubating for additional 1–2 min.

5. The resulting cell lysate can now be applied directly to a BD TALON CellThru Column, or the lysate supernatant can be applied to any other BD TALON Resin column after centrifuging at 10,000–12,000 \times g for 20 min at 4°C.

B. Standard Sample Preparation to Isolate Native Proteins

1. Harvest the cell culture by centrifugation at 1,000–3,000 \times g for 15 min at 4°C. Remove the supernatant. If yield is low, use the mild extraction method described in Step 6, below.
2. Resuspend the cell pellet by vortexing in 2 ml of chilled 1X Equilibration/Wash Buffer (4°C) per 25 ml of culture \leq 100 ml. For cultures >1 L, resuspend the pellet in 1–2% of the original culture volume.

Note: You may omit Steps 3–4 if lysozyme treatment interferes with your protein's functionality.

VI. Sample Preparation *continued*

3. Add lysozyme to the 1X Equilibration/Wash Buffer for a concentration of 0.75 mg/ml. To reduce the chance of introducing proteases, use the highest purity lysozyme available.

4. Incubate at room temperature for 20–30 min.

Note: Incubations at room temperature result in elevated proteolytic activities. Alternatively, you can use lysozyme at 4°C with lower efficiency. If this treatment hydrolyzes the target protein, use the method described in Step 6 (below). Alternatively, disrupt the cells by repeated freeze/thaw cycles; that is, flash-freezing the cell suspension in a dry ice-ethanol bath and thawing in chilled H₂O.

5. If your sample is ≤50 ml, sonicate it 3 x 10 sec, with a 30-sec pause on ice between each burst. If your sample is ≥ 200 ml, sonicate it 3 x 30 sec, with a 2 min pause on ice between each burst. Proceed to Step 7.

Note: Excessive sonication can destroy protein functionality.

6. **[Optional]: High-yield, mild extraction method.** Transfer the cells to a chilled mortar and grind 1 part cells with 2.5 parts Alumina (Sigma Cat. No.A-2039) for 2–3 min or until the composition of the mixture becomes paste-like. Add 2 ml chilled 1X Equilibration/Wash Buffer (4°C) per 25 ml culture.

Note: If there is a high level of proteolytic activity in the cell lysate, we recommend adding 1 mM EDTA (final concentration) to the Equilibration/Wash Buffer in order to inhibit metalloproteases during the extraction. Before application of the sample to the BD TALON Resin, EDTA must be removed by gel filtration chromatography (PD-10, Amersham, Pharmacia) equilibrated with the Equilibration/Wash Buffer for IMAC.

7. Centrifuge the cell extract at 10,000–12,000 x g for 20 min at 4°C to pellet any insoluble material.

8. Carefully transfer the supernatant to a clean tube without disturbing the pellet. This is the clarified sample.

9. Reserve a small portion of the clarified sample at 4°C for SDS/PAGE analysis.

10. If this is the first time you have prepared clarified samples from cells expressing a particular recombinant protein, we recommend that you estimate the protein's expression level in that host strain. To do so, perform a small-scale purification, and then analyze a portion by SDS/PAGE in parallel with protein standards. Once expression is observed, proceed with the appropriate purification protocol, below.

C. Standard Sample Preparation to Isolate Denatured Proteins

1. Harvest 20–25 ml of cell culture by centrifugation at 1,000–3,000 x g for 15 min at 4°C.

2. Resuspend the pellet in 2 ml of **Denaturing** 1X Equilibration/Wash Buffer (pH 7.0) per 20–25 ml of culture.

3. Gently agitate or stir the sample until it becomes translucent.

VI. Sample Preparation *continued*

4. Centrifuge the sample at 10,000–12,000 x g for 20 min at 4°C to pellet any insoluble material.
5. Carefully transfer the supernatant to a clean tube without disturbing the pellet. This is the clarified sample.
6. Set aside a small portion of the clarified sample for SDS/PAGE analysis. Then proceed with the appropriate purification protocol, below.

Note: Samples containing 6 M guanidinium must be dialyzed overnight against buffer containing 8 M urea before loading on a gel.

D. Standard Sample Preparation for BD TALON™ CellThru Resin

Native Protein

1. Harvest the cell culture by centrifugation at 1,000–3,000 x g for 15 min at 4°C. Remove the supernatant.
2. Resuspend the cell pellet by vortexing in 2 ml of chilled 1X Equilibration/Wash Buffer (4°C) per 25 ml of culture ≤100 ml. For cultures >1 L, resuspend the pellet in 1–2% of the original culture volume.

Note: You may omit Steps 3–4 if lysozyme treatment interferes with your protein's function.

3. Add lysozyme to the 1X Equilibration/Wash Buffer for a concentration of 0.75 mg/ml. To reduce the chance of introducing proteases, use the highest purity lysozyme available.
4. Incubate at room temperature for 20–30 min.

Note: Incubations at room temperature result in elevated proteolytic activities. Alternatively, you can use lysozyme at 4°C with lower efficiency. If this treatment hydrolyzes the target protein, use the method described in Step 6. Alternatively, disrupt the cells by repeated freeze/thaw cycles; that is, flash-freezing the cell suspension in a dry ice-ethanol bath and thawing in chilled H₂O.

5. If your sample is ≤50 ml, sonicate it 3 x 10 sec, with a pause for 30 sec on ice between each burst. If your sample is ≥ 200 ml, sonicate it 3 x 30 sec, with a 2 min pause on ice between each burst.

Note: Excessive sonication can destroy protein functionality.

6. Store a small portion of the clarified sample at 4°C for SDS/PAGE analysis.

Denatured Protein

1. Harvest 20–25 ml of cell culture by centrifugation at 1,000–3,000 x g for 15 min at 4°C.
2. Resuspend the pellet in 2 ml of **Denaturing** 1X Equilibration/Wash Buffer (pH 7.0) per 20–25 ml of culture.
3. Gently agitate or stir the sample until it becomes translucent.

VI. Sample Preparation *continued*

4. Set aside a small portion of the clarified sample for SDS/PAGE analysis. Then proceed with the appropriate purification protocol, below.

Note: Samples containing 6 M guanidinium must be dialyzed overnight against buffer containing 8 M urea before loading on a gel.

E. Standard HT 96-Well Sample Preparation

1. Grow cells in appropriate format for high throughput analysis.
2. Centrifuge if necessary and remove supernatant
3. If the target proteins are secreted in the medium, utilize the supernatant as a starting material and proceed to Step 6. If the target proteins are intracellular, proceed to the next step.
4. Disrupt the cells in presence of 1X Equilibration/Wash Buffer (Use 2 ml of buffer per 200 mg of cells per purification well).
5. Centrifuge extracts and collect the supernatant to be used as a starting material.
6. Remove 50 μ l of each sample for protein concentration analyses.

VII. Protein Purification Protocols

PLEASE READ ENTIRE PROTOCOL BEFORE STARTING.

A. General Information

1. Perform all manipulations at 4–8°C in order to maintain protein stability and improve yield.
2. **This protocol is designed using the Equilibration/Wash Buffer (pH 7.0).** If your target protein is more stable at pH 8.0, or if it does not adsorb at pH 7.0, use the Equilibration Buffer at pH 8.0 (instead of the Equilibration/Wash Buffer) during extraction and wash steps.
3. A reducing agent, such as 10 mM β -ME, or a protease inhibitor, such as PMSF, in the Equilibration/Wash Buffer (pH 7.0), may improve the structural stability of fragile proteins during sample preparation. See Appendix A for compatibility information.

Note: Depending on the concentration and volume of the additive you wish to use, you may need to remake the buffers to preserve the recommended concentration of NaCl and buffering agent. **DTT and DTE are not compatible with this BD TALON protocol in any concentration.**

4. If the cell lysate contains a high level of proteolytic activity, we recommend adding 1 mM EDTA to the Equilibration/Wash Buffer (pH 7.0) to inhibit metalloproteases during the extraction. **However, before applying the sample to the resin, remove EDTA using a gel**

VII. Protein Purification Protocols *continued*

filtration column (such as PD-10, Amersham, Pharmacia) equilibrated with the Equilibration/Wash Buffer. In some cases, the host cell produces low molecular weight chelators that can also be removed using gel filtration.

Chelators can be detected easily by applying your sample to a small column packed with BD TALON™ Resin. If the top of the column loses its characteristic pink color, and the colorless front moves in the direction of the flow, or if you obtain pink fractions during batch adsorption, you must equilibrate the sample using a gel filtration column.

- Overexpressed recombinant proteins can accumulate in insoluble inclusion bodies. In order to determine optimal extraction/purification conditions, you must determine the distribution of the protein in soluble and insoluble forms. Perform a preliminary SDS/PAGE analysis of protein extracts obtained under native conditions, followed by extraction of the residual proteins under denaturing conditions. Take care to use the same extraction volumes for both native and denaturing extracts, and run the cell extract before induction as a control in one lane to identify the target protein. Use of denaturing conditions is recommended only if the biological activity of the target protein would not be affected by denaturation or is unimportant. It is preferable to use native conditions for extraction even if only 5–10% of the target protein is soluble.
- The buffer volumes in the following protocols were optimized for purifying the HAT-DHFR fusion protein from 20–25 ml of *E. coli* culture. Depending on the expression level and anticipated yield, you may need to adjust the buffer volumes for other proteins. As a starting point, use 2 ml of buffer per 20–25 ml of culture.
- If you are purifying protein from harvested eukaryotic cells, lyse the cells in an appropriate buffer containing a mild detergent (Sambrook *et al.*, 1989). See Appendix A for compatible buffer additives. *Note that EDTA and EGTA are not compatible with these protocols because these metal-chelating reagents strip the cobalt from the resin.*
- Carefully check the sample's appearance after lysis or sonication. Bacterial samples often remain viscous from incomplete shearing of genomic DNA. Complete DNA fragmentation improves protein yields and allows efficient removal of cellular debris during centrifugation. You may decrease the sample's viscosity by digestion for 20–30 min at room temperature with 2.5 µg/ml of DNase I. Remember that proteolytic activity is much higher at room temperature. Alternatively, dilute the sample five-fold with Equilibration/Wash Buffer before applying it to the resin. This procedure should not significantly affect recovery.

VII. Protein Purification Protocols *continued*

Notes on Protein Purification methods using BD TALON™ Resin

The following general guidelines are used for purifying polyhistidine-tagged protein from transformed *E. coli* cultures. Table I provides an overview of BD TALON™ Resin protein purification methods and applications. Choose a method that best suits your research needs.

- Use 2 ml of resin suspension per ~3 mg of anticipated polyhistidine-tagged protein. 2 ml of homogeneously resuspended resin will provide 1 ml (bed volume) of BD TALON Resin.
- The buffers and purification conditions should work well for most soluble, monomeric proteins expressed in *E. coli*.
- Initially, test each different expression system and polyhistidine-tagged protein in small-scale batch purification to determine expression levels and to optimize the protocol. A mini-scale batch purification protocol is provided in Appendix B; alternatively, you can use a BD TALONspin column.
- Purification methods that work with Nickel or Zinc-based IMAC resins should also work with these resins. However, some optimization may be required.

BD TALON™ CellThru Considerations

The procedure for purifying polyhistidine-tagged proteins using BD TALON CellThru Resin is essentially the same as other BD TALON Resins with the following significant differences.

1. Extracellular Proteins

If there are no chelating agents in the fermentation liquid and the pH is ≥ 7.0 , you can apply sample directly onto a prepacked column. Otherwise, a desalting/equilibration step is necessary (such as ultrafiltration or gel filtration with Sephadex G25).

2. Intracellular Proteins

For purifying intracellular proteins, apply the sonicated sample containing your target proteins, directly onto a prepacked column. There is no need for centrifugation. Electrophoresis will reveal that some of the target protein has passed through the column without adsorption. To a large extent the material passing through the column is insoluble protein, which would normally have been removed during high-speed centrifugation. The amount of nonadsorbed target protein will also vary as a function of sonication efficiency.

3. Chromatography Considerations

BD TALON CellThru Beads have a diameter of 300–500 μm ; therefore, use a column with a filter pore size of 90–130 μm to adequately pass cellular debris. We recommend using our CellThru 2-ml & 10-ml Disposable Columns (Cat. No.8914-1 & 8915-1). The 2-ml columns are

VII. Protein Purification Protocols *continued*

suitable for 1–2 ml bed volumes, while the 10-ml columns are suitable for 5–10 ml bed volumes. Because the column filters have a larger pore size and permit higher flow rates, you may need to incubate your sample with the resin for 5 minutes before letting it flow through. If necessary, pass the sample through the column a second time.

The technique of expanded bed chromatography works well with these resins as the material can flow through the resin more effectively. Flow rates may have to be adjusted to get the maximum binding efficiency when using this technique.

B. Batch/Gravity-Flow Column Purification

For column IMAC using BD TALON Resins, we recommend a hybrid batch/gravity-flow procedure. This method combines the speed and convenience of a batch procedure with the exceptionally high purity of the gravity-flow column method. In this hybrid procedure, the binding and initial washing steps are performed in a batch format to save time, eliminate extraneous debris, and avoid column clogging. After the initial washes, the resin is transferred to a column for additional washing and protein elution.

1. Thoroughly resuspend the BD TALON Resin.
2. Immediately transfer the required amount of resin suspension to a sterile tube that will accommodate 10–20 times the resin bed volume.
3. Centrifuge at 700 x g for 2 min to pellet the resin.
4. Remove and discard the supernatant.
5. Add 10 bed volumes of 1X Equilibration/Wash Buffer and mix briefly to pre-equilibrate the resin.
6. Re-centrifuge at 700 x g for 2 min to pellet the resin. Discard the supernatant.
7. Repeat Steps 5 and 6.
8. Add the clarified sample from Section VI.A, B or C to the resin.
9. Gently agitate at room temperature for 20 min on a platform shaker to allow the polyhistidine-tagged protein to bind the resin.
10. Centrifuge at 700 x g for 5 min.
11. Carefully remove as much supernatant as possible without disturbing the resin pellet.
12. Wash the resin by adding 10–20 bed volumes of 1X Equilibration/Wash Buffer. Gently agitate the suspension at room temperature for 10 min on a platform shaker to promote thorough washing.
13. Centrifuge at 700 x g for 5 min.
14. Remove and discard the supernatant.
15. Repeat Steps 12–14.

VII. Protein Purification Protocols *continued*

16. Add one bed volume of the 1X Equilibration/Wash Buffer to the resin, and resuspend by vortexing.
17. Transfer the resin to a 2-ml gravity-flow column with an end-cap in place, and allow the resin to settle out of suspension.
18. Remove the end-cap, and allow the buffer to drain until it reaches the top of the resin bed, making sure no air bubbles are trapped in the resin bed.
19. Wash column once with 5 bed volumes of 1X Equilibration/Wash Buffer.
20. **[Optional]:** If necessary, repeat Step 19 under more stringent conditions using 5–10 mM imidazole in 1X Equilibration/Wash Buffer (Section IV.D).
21. Elute the polyhistidine-tagged protein by adding 5 bed volumes of Elution Buffer to the column. Collect the eluate in 500- μ l fractions.
Note: Under most conditions, the majority of the polyhistidine-tagged protein will be recovered in the first two bed volumes.
22. Use spectrophotometric and SDS/PAGE analyses to determine which fraction(s) contain(s) the majority of the polyhistidine-tagged protein.
Note: Use a Bradford protein assay (Bradford, 1976) or UV absorbance at 280 nm. Use UV absorbance only if you are eluting sufficient protein to exceed the absorbance of the imidazole at 280 nm. Alternatively, dialyze the fractions overnight against the Equilibration/Wash Buffer, and then measure their UV absorbance at 280 nm.

C. Large-Scale Batch Purification

This method purifies polyhistidine-tagged proteins faster than gravity-flow columns; however, batch washes remove impurities less efficiently than gravity-flow columns. Therefore, they require larger wash buffer volumes to obtain pure polyhistidine-tagged proteins.

1. Thoroughly resuspend BD TALON Resin.
2. Transfer required amount of resin to a glass filter with a pore size of 10–20 μ m.
3. Apply a vacuum to the filter to remove excess ethanol.
4. Add 5 bed volumes of deionized water to the resin, and apply vacuum.
5. Add 5 bed volumes of 1X Equilibration/Wash Buffer to the resin, and apply vacuum.
6. Repeat Step 5 two times.
7. Add crude lysate (CellThru Resin) or clarified sample (other than CellThru Resin) to the resin, and mix for 3–5 min.
8. Apply vacuum and collect the filtrate.
9. Wash the resin by adding 10–20 bed volumes of 1X Equilibration/Wash Buffer. Gently agitate the suspension at room temperature for 10 min on a platform shaker to promote thorough washing.
10. Apply vacuum to remove buffer.

VII. Protein Purification Protocols *continued*

11. Repeat the above wash (Steps 9–10) 2–3 times.
12. **[Optional]:** If necessary, repeat Step 11 under more stringent conditions using 5 mM imidazole in 1X Equilibration/Wash Buffer (Section IV.D).
13. Elute the polyhistidine-tagged protein by adding 5 bed volumes of Elution Buffer.
14. Gently agitate suspension at room temperature for 5 min.
15. Apply vacuum, and collect the purified polyhistidine-tagged protein.
16. Repeat Steps 13–15 two times, collecting separate fractions.
17. Use spectrophotometric and SDS/PAGE analyses to determine which fraction(s) contain(s) the majority of the polyhistidine-tagged protein.
Note: Samples containing 6 M guanidinium must be dialyzed overnight against buffer containing 8 M urea before loading on a gel.

D. Medium-Pressure Column Purification

1. Assemble column according to the manufacturer's instructions.
2. Thoroughly resuspend BD TALON Superflow Resin. Slowly pour the slurry into the column, and avoid introducing air bubbles.
3. Allow resin to settle. Accelerate this process by allowing the buffer to flow through the column with a peristaltic pump attached to the output of the column. Do not exceed a flow rate of 5 ml/min/cm². Do not allow the resin to dry out. If this occurs, resuspend the resin and repack the column.
4. Insert and adjust the top adaptor and connect the column to the chromatography system according to manufacturer's instructions.
Note: Avoid trapping air between the adaptor and the resin surface.
5. Equilibrate the column with 1X Equilibration/Wash Buffer. Do not exceed a 5 ml/min/cm² flow rate. Monitor the eluant at 280 nm; the baseline should be stable after washing with 5–10 column volumes.
6. Apply the clarified sample to the column after filtering it through a 0.22- μ m filter and wash with Equilibration/Wash Buffer until the baseline (280 nm) is stable. Monitor column backpressure during sample application. Start collecting fractions.

Note: If the sample is very viscous, the column pressure may exceed the recommended value (150 psi, 1.0 MPa). Reduce the flow rate or dilute the sample to bring the pressure into an acceptable range.

Load the sample at a flow rate of 0.5–1.0 ml/min/cm² to ensure that the polyhistidine-tagged protein will bind to the resin. If the protein does not bind, reduce the flow rate further. If desired, increase the flow rate for washing and protein elution.

If the target protein is unstable at room temperature, perform the chromatography at 4°C. Alternatively, use flow rates up to 5 ml/min/cm²

VII. Protein Purification Protocols *continued*

to load, wash, and elute the protein. Capacity will decrease by 10–15%, but on average, a chromatography run should only take 15–20 min.

7. Wash column with 10–20 column-volumes of Equilibration/Wash Buffer, or until the baseline at 280 nm is stable. If necessary, add 5–10 mM imidazole to the Equilibration/Wash Buffer.
8. Elute the polyhistidine-tagged protein with 5–10 column-volumes of Elution Buffer. The polyhistidine-tagged protein usually elutes in the second and third column volumes.
9. Use spectrophotometric and SDS/PAGE analyses to determine which fraction(s) contain(s) the majority of the polyhistidine-tagged protein.
Note: Samples containing 6 M guanidinium must be dialyzed overnight against buffer containing 8 M urea before loading on a gel.
10. If you plan to store, regenerate, and reuse a resin-packed column, see Section VIII.C.

E. BD TALONspin™ Column Purification

Important Points

- Before proceeding with purification, determine the concentration of polyhistidine-tagged protein in your sample using the mini-batch screening protocol (Appendix B). Alternatively, run a sample of the clarified lysate directly on SDS/PAGE, and estimate the amount of polyhistidine-tagged protein by band intensity.
 - Avoid excessively concentrated or viscous lysates. See Troubleshooting (Section IX.B.2) for tips on reducing sample viscosity.
 - If the concentration of polyhistidine-tagged protein in the lysate is very dilute, use one column to enrich the protein from several 0.6–1 ml lysate aliquots. Simply repeat Steps 11–16 (below) until the desired amount of lysate has been processed. Alternatively, concentrate the polyhistidine-tagged protein by reducing the sample volume.
 - The centrifugation rotor and speed may affect your results. Ideally, you should centrifuge BD TALONspin Columns in a swinging bucket rotor to allow the sample to pass through the resin uniformly. However, a fixed angle rotor or a microcentrifuge is also acceptable. Centrifugation speeds higher than 700 x g may cause irregularities in the flow of solution through the resin bed, and thus, decrease the performance of the column.
1. Hold the BD TALONspin Column upright and flick it until all resin falls to the bottom of the column.
 2. Snap off the breakaway seal.
 3. Place column in the 2-ml microcentrifuge tube.
 4. Save white end-cap for later use.
 5. Remove the clear top-cap and centrifuge column at 700 x g for 2 min

VII. Protein Purification Protocols *continued*

to remove the storage buffer from the resin bed.

Note: The resin bed will appear semi-dry after centrifugation.

6. Remove column from centrifuge, and place the white end-cap over the male luer fitting.
7. Add 1 ml 1X Equilibration/Wash Buffer and mix briefly to pre-equilibrate the resin.
8. Re-centrifuge at 700 x g for 2 min to pellet the resin. Discard the supernatant.
9. Repeat Steps 7 and 8, twice.
10. Add the clarified sample from Section VI.A, B or C to the resin.
11. Add 0.6–1 ml of sample to the column, and replace the clear top-cap.
12. Allow sample to passively wet the resin bed for 30 sec.
13. Mix or vortex contents briskly for 1–2 sec, completely resuspending the resin in the lysate.
14. Gently agitate the suspension for 5 min to allow polyhistidine-tagged protein binding. **Do not vortex.**
15. Remove both caps from column and place column inside the 2-ml microcentrifuge tube.
16. Centrifuge at 700 x g for 2 min.
17. Remove the column and microcentrifuge tube from the centrifuge rotor, making sure that all of the sample has passed through the resin bed.
Note: Viscous samples may require additional centrifugation.
18. Save the 2-ml tube, but discard the flowthrough.
19. Place microcentrifuge tube in rotor.
20. Place white end-cap on the column, and add 1 ml of 1X Equilibration/Wash Buffer. Close the column with the clear top-cap.
21. Allow the buffer to passively wet the resin bed for 30 sec.
22. Agitate or vortex briskly for a few seconds until the resin is completely resuspended.
23. Gently agitate for 5 min.
24. Remove both caps, and centrifuge at 700 x g for 2 min.
25. Repeat Steps 18–24. Repeat twice for particularly concentrated lysates, or if necessary, to improve purity.
26. Examine the resin bed to ensure that it appears semi-dry, and to ensure that all wash buffer has drained from the resin bed and the column end.
27. Discard the used 2-ml microcentrifuge tube.
28. If necessary, repeat the spin to remove all traces of wash buffer.
29. Replace the white end-cap on the spin column.

VII. Protein Purification Protocols *continued*

30. Add 400–600 μ l of Elution Buffer.

Note: Alternatively, use 100 mM EDTA (pH 8.0) if it does not interfere with downstream applications of the protein. Samples eluted with EDTA will also contain cobalt.

31. Allow 1 min for Elution Buffer to passively wet the resin bed.

32. Briefly agitate or vortex to resuspend the resin.

33. Place a fresh 2-ml collection tube into centrifuge rotor.

34. Remove both caps and place column into the 2-ml collection tube.

35. Centrifuge sample at 700 x g for 2 min.

36. Repeat Steps 30–35.

Note: The polyhistidine-tagged protein sample can generally be recovered in 800–1,200 μ l of Elution Buffer, but it may be necessary to use a larger Elution Buffer volume or repeat Steps 30–35.

37. Determine polyhistidine-tagged protein yield using gel or spectrophotometric analysis.

Note: If the purity of the polyhistidine-tagged protein preparation is unsatisfactory, refer to the procedure in the Troubleshooting Guide Section IX.C.2.

F. BD TALON™ HT 96-Well Purification Protocol

Each well of the BD TALON HT 96-Well Plate has a capacity of up to 1.0 mg of polyhistidine tagged protein. In order to obtain the maximum yield of pure protein, do not attempt to load more than 1.0 mg of polyhistidine-tagged protein/well. Also, **observe the following guidelines:**

- When using pipette tips to mix the resin, use wide-bore tips, or cut the tips to make the opening wider. This will reduce mechanical damage to proteins as well as resin.
- BD TALON Resin is designed to permit buffer to flow through freely. Therefore, when the HT 96-Well Plate is not on the vacuum manifold or over a Collection Plate, we recommend that it is kept on the Base Seal that acts as a temporary stopper.
- The amount of sample applied to a well should not exceed the capacity of 1.0 mg/well.
- Avoid overdrying the resin under the vacuum. For the best results, keep the resin wet.
- When using vacuum manifold, adjust the vacuum to obtain a flow rate of 1–2 drops per sec (~100–200 mm Hg or 2–4 psi).

1. Unpacking and removal of seals

HT 96-Well plates come with solid plate seals to prevent resin from leaking during transportation. Before removing the upper seal, we recommend performing a 2 min centrifugation step at 500 x g to pack resin that might have adhered to the silicon lid during transportation.

VII. Protein Purification Protocols *continued*

After this procedure, the upper seal can be removed and the steps outlined in the purification protocols can be performed.

If you do not desire to use all 96 wells, the plate seals can be cut so that only the wells that are needed are exposed. After chromatography, the removed portion of the plate seals can be replaced and the plate can be stored at 4°C until the remaining wells are used. When stored, the resin in unused wells should be covered in 20% ethanol.

2. HT 96-Well Plate Equilibration

- a. Remove the top and bottom seals from the plate.
- b. Place the plate on the manifold and apply vacuum to remove storage solution or centrifuge 5 min at 700 X g.
- c. Add 1 ml of deionized water to each well of the plate. Apply vacuum or centrifuge to drain the water from the wells. Repeat twice.
- d. Add 1 ml of 1X Equilibration/Wash Buffer to each well of the plate and apply vacuum or centrifuge to drain the buffer from the well. Repeat twice.

3. Vacuum Purification

When performing vacuum purification, adjust the vacuum to obtain a flow rate of 1–2 drops per sec (~100–200 mm Hg or 2–4 psi). In addition, avoid overdrying the resin which introduces air bubbles and reduces performance.

- a. Apply 1.5 ml of the starting sample (See Section VI.E) per well. Mix the sample with the resin shortly by vortexing the plate or pipetting up and down inside the wells. Leave the plate on ice for 5–10 min mixing samples every 2 min.
- b. Place the plate on the vacuum manifold, apply vacuum and let the excess liquid drain into the manifold. Firmly press all four sides of the plate to the rubber gasket of the vacuum manifold. Ensure by observation that all wells have been drained of buffer.
- c. Repeat Steps 3.a and 3.b if additional loading is necessary.
- d. Add 1 ml of 1X Equilibration/Wash Buffer and suspend the resin by vortexing the plate or pipetting up and down inside the wells.
- e. Place the plate on the vacuum manifold, apply vacuum and let the excess liquid drain into the manifold. Firmly press all four sides of the plate to the rubber gasket of the vacuum manifold. Ensure by observation that all wells have been drained of buffer.
- f. Repeat Steps 3.d and 3.e twice.
- g. Add 1 ml of HT 96-Well Plate Wash Buffer (See Section III) to each well and suspend the resin by vortexing the plate or pipetting up and down inside the wells.

VII. Protein Purification Protocols *continued*

- h. Place the plate on the vacuum manifold, apply vacuum and let the excess liquid drain into the manifold. Firmly press all four sides of the plate to the rubber gasket of the vacuum manifold. Ensure by observation that all wells have been drained of buffer.
- i. Repeat Steps 3.g and 3.h five times.
- j. Remove the HT 96-Well Plate from the vacuum manifold. Drain the collected filtrate from the vacuum manifold.
- k. Place a Collection Plate inside the vacuum manifold and place the HT 96-Well Plate on the vacuum manifold.
Note: Before Eluting, ensure that the plate is over a Collection Plate or on the base seal.
- l. Add 200 μ l of 1X Elution buffer (Section III.A) and suspend the resin by vortexing the plate or pipetting up and down inside the wells.
- m. Place the HT 96-Well Plate on the vacuum manifold, apply vacuum, and let the eluate drain into the Collection Plate.
- n. Repeat elution (Steps 3.m and 3.n) twice.
- o. Determine amount of loaded and adsorbed protein in each well by Bradford Assay (Bradford, 1976).

4. Centrifuge Purification

Due to the variety of rotors and centrifuges that can be used, the following instructions are only general guidelines for successful purification:

- Do not utilize centrifugal force higher than 700 x g.
 - Ensure proper balance of the HT 96-Well Plate/Collection Plate inside the rotor.
 - When performing the centrifuge procedure below, extra Collection Plates are recommended. See Additional Materials Required for information on obtaining compatible plates.
- a. Add 1.5 ml of the starting sample per well (See Section VI.E). Mix the sample with the resin briefly by vortexing the plate or pipetting up and down inside the wells. Leave the plate on ice for 5–10 min mixing samples every 2 min.
 - b. Centrifuge the plate for 5 min. Ensure by observation that all wells have been drained of buffer.
 - c. Repeat Steps 4.a and 4.b if additional loading is necessary.
 - d. Add 1 ml of 1X Equilibration/Wash Buffer and suspend the resin by vortexing the plate or pipetting up and down inside the wells.
 - e. Centrifuge the plate for 5 min. Ensure by observation that all wells have been drained of buffer.
 - f. Repeat Steps 4.d and 4.e twice.

VII. Protein Purification Protocols *continued*

- g. Add 1 ml of HT 96-Well Plate Wash buffer and suspend the resin by vortexing the plate or pipetting up and down inside the wells.
- h. Centrifuge the HT 96-Well Plate for 5 min. Ensure by observation that all wells have been drained of buffer.
- i. Repeat Steps 4.g and 4.h five times.
- j. Drain the collected filtrate from the Collection Plate or use a fresh Collection Plate (See Additional Materials).
- k. Place the HT 96-Well Plate on the Collection Plate in the rotor and centrifuge 5 min. Ensure by observation that all wells have been drained of buffer.
Note: Before Eluting, ensure that the plate is over a Collection Plate or on the base seal.
- l. Add 200 μ l of 1X Elution buffer (Section III.A) and suspend the resin by vortexing the plate or pipetting up and down inside the wells.
- m. Centrifuge the HT 96-Well Plate on the Collection Plate for 5 min. Ensure by observation that all wells have been drained of buffer.
- n. Repeat elution (Steps 4.l and 4.m) twice.
- o. Determine amount of loaded and adsorbed protein in each well by Bradford Assay (Bradford, 1976).

VIII. Resin Washing, Reuse, and Storage

Generally, reuse BD TALON™ Resins 3–4 times before discarding or complete regeneration. The exact number of uses varies among preparations and application because of differences in redox potential, organic complexity, and debris content. To avoid possible cross-contamination, use a particular aliquot of resin to purify a single type of polyhistidine-tagged protein.

Important precautions

- BD TALONspin™ Columns are not reusable.
- Do not store BD TALON Resin in denaturants such as 6 M guanidinium.
- Do not store BD TALON Resin with bound imidazole: the resin should be washed with MES Buffer (pH 5.0) described in Section III: Additional Materials Required before reuse to remove the bound imidazole.

A. Stringent Wash (optional)

1. Wash resin with four bed volumes of 6 M guanidinium (pH 5.0) + 1% nonionic detergent.
2. Rinse resin with five bed volumes of distilled H₂O.
3. Store resin at 4°C in 20% nonbuffered ethanol containing 0.1% azide.

VIII. Resin Washing, Reuse, and Storage *continued*

B. Removing Imidazole

1. Wash resin with five bed volumes of 20 mM MES Buffer (pH 5.0) containing 0.1 M NaCl.
2. Rinse resin with five bed volumes of distilled H₂O.
3. Store resin at 4°C in 20% nonbuffered ethanol containing 0.1% azide.

C. Regeneration of Superflow Columns

Purification of polyhistidine-tagged proteins using imidazole gradients will cause the column to take on a purplish hue. Washing the column with 5–10 column volumes of 20 mM MES Buffer (pH 5.0) will restore the normal pink color and bring the absorbance at 280 nm back to the original baseline level. After equilibrating the column with Equilibration/Wash Buffer, the column is ready for reuse.

D. Complete Regeneration

1. Strip the resin of cobalt ions by washing with 10 bed volumes of 0.2 M EDTA (pH 7.0).
2. Wash excess EDTA from the resin with an additional 10 bed volumes of double distilled H₂O (ddH₂O).
3. Charge the resin with 50 mM CoCl₂ solution (10 bed volumes).
4. Wash resin with 7 bed volumes of ddH₂O followed by 3 bed volumes of 300 mM NaCl, and 3 bed volumes of ddH₂O to remove excess cobalt metal ions.
5. Equilibrate the resin with equilibration/wash buffer (10 bed volumes).
6. If you plan to use β-mercaptoethanol in subsequent buffers/procedures, then re-equilibrate the resin as follows before proceeding with further purifications:
 - a. Wash the resin with at least two bed volumes of Equilibration/Wash Buffer.
 - b. Re-equilibrate the resin with Equilibration/Wash Buffer containing β-mercaptoethanol.

IX. Troubleshooting Guide

A. Protein Expression

1. No expression

- **Bad vector construct**
Check sequence of the vector.
- **Bad transformation**
Make a plasmid miniprep and confirm sequence.
- **No inducing agent added to culture to induce expression**

2. Apparent low expression

- **Insoluble over-expressed protein**
Use denaturing extraction and purification conditions or reduce expression levels by lowering the amount of inducer.
- **Unsuitable expression conditions**
Check cell growth and inducer concentration; check for wild-type (nontransformed) or antibiotic resistant cells.
- **Protein is secreted**
Use fermentation liquid as starting sample for IMAC after proper buffering.

B. Loading/Washing

1. Polyhistidine-tagged protein elutes in the wash buffer

- **Problems with vector construction**
Ensure that protein and tag are in frame.
- **Buffer is not optimal**
Check the pH and composition of all buffers. Use a lower stringency wash buffer for all washing steps. For example, slightly increase the pH of the wash buffer or lower its imidazole concentration.
- **Protein degraded during extraction**
 - a. Use mild extraction conditions in the presence of protease inhibitors (e.g., β -ME and EDTA) at 4°C. **Be sure to remove EDTA before applying to BD TALON™ Resin.**
 - b. Make C-terminal construct.
 - c. Work quickly at 4°C to reduce the time for initial purification steps.
- **Reagent interferes with binding**
See Appendix A for reagent compatibilities. Dilute an aliquot of lysate (1:10), or sonicate, and check binding on a small scale. Try using a different polyhistidine-tagged protein as a control.

IX. Troubleshooting Guide *continued*

- **Tag is not accessible under native conditions**

If the protein fails to bind under native conditions, treat a small aliquot (<1 ml) with 6 M guanidinium and bind to 50 µl of resin. Then follow the mini-scale procedure in Appendix B. If the target protein binds to the resin under the denaturing conditions, then try to move the tag to the other terminus of the protein where it may be more exposed.

2. High back pressure during load of sample

- **High viscosity due to presence of DNA**

Use DNase I or dilute sample five fold (Section VII.A.8).

C. Elution

1. High amount of co-eluted impurities

- **Insufficient wash**

Use larger volumes of Equilibration/Wash Buffer.

- **Buffer compositions are not optimal**

- a. Check buffers used for sample preparation and wash steps.
- b. Check pH. The Equilibration/Wash Buffer should be pH 7.0. Contaminants will co-elute in buffers <pH 7.0.
- c. Increase volume of wash buffer and continue to wash resin bed until the A_{280} drops to zero.
- d. Increase counterion concentration up to 0.5 M NaCl or KCl to inhibit nonspecific ionic interactions.
- e. Add ethylene glycol or glycerol to inhibit nonspecific hydrophobic interactions.
- f. Add small amounts of nonionic detergent(s); this is particularly important when isolating proteins from a eukaryotic expression system.
- g. Add 5–10 mM imidazole to the Equilibration/Wash Buffer and use it as an intermediate wash step before elution.

- **Proteolytic product**

Use mild extraction conditions in presence of protease inhibitors (e.g., β -ME and EDTA) at 4°C. **Remove EDTA before applying to BD TALON™ Resin.**

- **Covalent attachment (Cys-Cys) of impurities to the protein**

Use 5–10 mM of β -ME in the Equilibration/Wash of Buffer.

- **Co-purifying histidine rich proteins**

- a) For HAT proteins, use enterokinase to remove HAT tag and rerun IMAC with mixture. Target protein will pass through the column, while impurities and tag will be adsorbed.

IX. Troubleshooting Guide *continued*

Note: Remove chelating ligands by gel filtration before loading the proteolytic mixture onto BD TALON Resin.

- b) Use second purification principle, such as size exclusion, ion exchange, hydrophobic, or thiophilic chromatography.

- **Protein sample is too concentrated and/or viscous**

Dilute sample 1:5 or 1:10 with additional buffer and centrifuge again before proceeding. Also, see the note on reducing sample viscosity after sonication in Section VII.A.8.

2. Excessive background after BD TALONspin™ Column procedure

- **Sample is too viscous**

- a) Reduce the viscosity of the sample (Section VII.A.8).
- b) Dilute clarified sample with an equal volume of Equilibration/Wash Buffer and load as two aliquots.
- c) Increase number of 1-ml washes.
- d) Use Equilibration/Wash Buffer (pH 7.0).
- e) Add 1–5 mM imidazole to Equilibration Buffer, pH 8.0 and use it as an intermediate wash step before elution.
- f) To re-purify a sample, perform the following after Step 37 in Section VII.E:
 1. Add 4 volumes of Equilibration/Wash Buffer to semi-purified sample.
 2. Load sample onto another BD TALONspin™ Column.
 3. Wash twice with 1 ml of Equilibration/Wash Buffer.
 4. Elute as before (Steps VIIE.30–35).

3. Column ceases to flow

- **Filter is clogged with subcellular debris**

Change column filters and centrifuge sample at 12,000 x g for 20–30 min at 4°C.

- **Proteins precipitated on the column**

Use a mild detergent such as Decanoyl-N-methylglucamide (MEGA-10, Sigma Cat. No.D6277) in the Equilibration/Wash Buffer.

- **The lower resin bed support may be clogged with cellular debris**

- a) Remove resin from clogged column and resuspend. Then wash it in a batch format and transfer to a fresh column.
- b) Use a syringe filled with wash buffer or reverse the pump on the column to gently run the column backwards. In addition, test for tubing blockages in a similar manner. Apply gentle pressure. Do not exceed a 1 drop/sec flow rate.

IX. Troubleshooting Guide *continued*

4. Polyhistidine-tagged proteins do not elute

- **Elution Buffer is less than optimal**

- a) Elute with 150 mM imidazole or pH 4.0 buffer.
- b) For really tough elution problems, you can strip off the protein using 100 mM EDTA (pH 8.0); however, doing so will remove the cobalt from the resin and deposit it in your protein sample.
- c) Add 1–5 mM β -ME to reduce disulfide linkages. Supplement buffer with 1% nonionic detergent.
- d) Purify polyhistidine-tagged protein under denaturing conditions.

D. Changes in Resin

1. Loss of Co^{2+}

- **Presence of chelators in sample**

- Remove chelators from sample by gel filtration
Regenerate resin as described in Section VIII.D.

2. Gray or brown resin

- **BDTALON™ Resin exposed to reducing agents or high concentration of β -ME**

- Completely remove reducing agents, such as DTE or DTT, or if possible, by gel filtration with β -ME. Reduce β -ME concentration (≤ 5 mM).

3. Resin particles aggregate or exhibit change in consistency

- **DNA cross-linking**

- a) Increase ionic strength of the buffers by using ≤ 500 mM NaCl or KCl.
- b) Vigorously sonicate sample to shear DNA.
- c) Pretreat sample with 100 $\mu\text{g/ml}$ DNase I at 30°C for 30 min.
- d) Dilute sample 1:5–1:10 with buffer, and repeat.
- e) Avoid long-term storage in denaturants.

E. Analysis

1. High background on silver-stained gels

- **Nucleic acid**

- a) Supplement buffer with 0.2–0.5 M NaCl or KCl. Repeat purification.
- b) Shear DNA more vigorously.
- c) Use DNase I in the extraction procedure.

IX. Troubleshooting Guide *continued*

2. Nonfunctional protein

- **Protein was damaged by sonication**

- a) Conduct a time-course assay to determine the minimum sonication time needed to disrupt the cells while maintaining the native protein/enzyme function. For example, sonicate samples at a medium-high setting for 0, 20, and 30 sec. Then perform protein or enzyme function assays and measure the A_{280} of each sample.
- b) Perform the lysis or sonication procedure on ice.

F. Reuse

1. Binding drops below original capacity

- **Lysate contains naturally occurring reducing agent or a non-specific polyanion may be obscuring the metal binding sites.**

- a) Use a larger volume of previously used resin.
- b) Replace used resin with fresh resin.
- c) Wash resin with 6 M guanidinium (pH 5.0) + 1% Triton X-100 or SDS, and re-equilibrate before use.

X. References

- Bradford, M. M. (1976) A rapid and sensitive method for the quantification of microgram quantities of protein utilizing the principle of protein-dye binding. *Anal. Biochem.* **72**:248–254.
- Bush, G. L., Tassin, A.-M., Friden, H. & Meyer, D. I. (1991) Secretion in yeast: purification and *in vitro* translocation of chemical amounts of prepro-alpha-factor. *J. Biol. Chem.* **266**:13811–13814.
- Chaga, G., Bochkariov, D. E., Jokhadze, G. G., Hopp, J. & Nelson, P. (1999) Natural poly-histidine affinity tag for purification of recombinant proteins on cobalt(II)-carboxymethylaspartate crosslinked agarose. *J. Chromatogr.* **864**:247–256.
- Hemdan, E. S. & Porath, J. (1985a). Development of Immobilized Metal Affinity Chromatography II. Interaction of amino acids with immobilized nickel iminodiacetate. *J. Chromatogr.* **323**:255–264.
- Hemdan, E. S. & Porath, J. (1985b). Development of Immobilized Metal Affinity Chromatography III. Interaction of oligopeptides with immobilized nickel iminodiacetate. *J. Chromatogr.* **323**:265–272.
- Hochuli, E., Döbeli, H. & Schacher, A. (1987) New metal chelate adsorbent selective for proteins and peptides containing neighboring histidine residues. *J. Chromatogr.* **411**:177–184.
- Hochuli, E., Bannwarth, W., Döbeli, H., Gentz, R. & Stüber, D. (1988) Genetic approach to facilitate purification of novel recombinant proteins with a novel metal chelate adsorbent. *Bio/Technology* **6**:1321–1325.
- Kasher, M. S., Wakulchik, M., Cook, J. A. & Smith, M. (1993) One-step purification of recombinant human papillomavirus Type 16 E7 oncoprotein and its binding to the retinoblastoma gene product. *BioTechniques* **14**:630–641.
- Porath, J. (1985). Immobilized Metal Ion Affinity Chromatography—A Powerful Method for Protein Purification. In *Modern Methods in Protein Chemistry*, (pp. 85–95). H. Tschelsche (Ed.), (Berlin & NY: Walter de Gruyter & Co).
- Porath, J., Carlsson, J., Olsson, I. & Belfrage, G. (1975). Metal chelate affinity chromatography, a new approach to protein fractionation. *Nature* **258**:598–599.
- Sambrook, J., Fritsch, E. F. & Maniatis, T. (1989) *Molecular Cloning: A Laboratory Manual* (Cold Spring Harbor Laboratory, Cold Spring Harbor, NY).
- Sulkowski, E. (1985). Purification of proteins by IMAC. *Trends Biotechnol.* **3**:1–7.
- Zhao, Y.-J., Sulkowski, E. & Porath, J. (1991). Surface topography of histidine residues in lysozymes. *Eur. J. Biochem.* **202**:1115–1119.

XI. Related Products

For a complete listing of all BD Biosciences Clontech products, please visit www.bdbiosciences.com/clontech

	Cat. No.	New Cat. No.
• CellThru 2-ml Disposable Columns	8914-1	635512
• CellThru 10-ml Disposable Columns	8915-1	635513
• 6xHis Monoclonal Antibody (Albumin Free)	8916-1	631212
• 6xHN Polyclonal Antibody	8940-1	631213
• BD PRO™Tet 6xHN Bacterial Expression System	K1628-1	631203
• HAT Protein Expression and Purification System	K6050-1	631205
• HAT Polyclonal Antibody	8909-1	631211
• pHAT20 Vector	8921-1	631202
• Glutathione-Uniflow Resin	8912-1	635610
	8912-2	635611
• Glutathione-Superflow Resin	8911-1	635607
	8911 -2	635608
• GST Purification Kit	K1251-1	635619
• Thiophilic Uniflow	8913-1	635613
	8913-2	635614
• Thiophilic Superflow	8917-1	635616
	8917-2	635617
• NucleoVac	4071-1	636030

Appendix A. Reagent Compatibilities and Incompatibilities

A. Compatible reagents

Table III shows the maximum concentrations of each reagent tested at BD Biosciences Clontech. Higher levels may be acceptable, but they should be tested before use. Note that some of these reagents may partially or completely denature your protein.

TABLE III. REAGENT COMPATIBILITY

Reagent	Acceptable Concentration
β -Mercaptoethanol ^a	10 mM (with caution)
CHAPS ^b	1% (with caution)
Ethanol ^c	30%
Ethylene glycol	30%
HEPES	50 mM
Glycerol	20%
Guanidinium ^a	6 M
Imidazole ^d	200 mM at pH 7.0–8.0, for elution
KCl	500 mM
MES	20 mM
MOPS	50 mM
NaCl	1.0 M
NP-40	1%
SDS ^b	1% with caution
TRIS ^e	50 mM
Urea	8 M

^a Use resin immediately after equilibrating with buffers containing these reagents. Otherwise, the resin will change color. Do not store resin in buffers containing these reagents.

^b Ionic detergents like CHAPS (3-[(30Cholamidopropyl)-dimethylammonio]-1-propane-sulfonate), SDS (sodium dodecyl sulfate), and Sarkosyl are compatible up to 1%. However, due to their charged nature, you should anticipate interference with binding.

^c Ethanol may precipitate proteins, causing low yields and column clogging.

^d Imidazole cannot be used at concentrations higher than 5–10 mM for loading polyhistidine-tagged proteins, because it competes with the histidine side chains (imidazole groups) for binding to the immobilized metal ions.

^e TRIS coordinates weakly with metal ions, causing a decrease in capacity.

B. Incompatible reagents

These reagents are incompatible at any concentration:

- DTT (dithiothreitol) and DTE (dithioerythritol)
- EDTA (ethylenediaminetetraacetic acid) and EGTA (ethylene glycol-bis([β -amino-ethyl ether])

Note: Although you can use EDTA at indicated points, it must be removed from the sample by gel filtration prior to applying it to BD TALON™ Resins.

Appendix B. Mini-Scale Protein Purification Protocol

Mini-scale protein purification is ideal for any of the following:

- (a) checking for a polyhistidine-tagged protein
- (b) determining expression levels
- (c) testing buffer conditions

You can also use a BD TALONspin™ Column with this procedure.

We recommend that you set aside a sample after each critical step of the procedure, and analyze all samples by SDS/PAGE.

Important

- This protocol is not intended for obtaining highly purified polyhistidine-tagged protein samples. Furthermore, protein samples eluted with EDTA (Step 19, below) will contain cobalt and EDTA, which may seriously inhibit enzyme activity and may cause the protein to precipitate.
- This protocol was optimized using denaturing conditions at pH 8.0. If you wish to obtain native samples, then substitute buffers accordingly. You may also need to use lysozyme (0.75 mg/ml of native buffer) to completely disrupt the cells in Step 5.

-
1. Transfer 1 ml of expression culture to a 1.5-ml microcentrifuge tube.
 2. Centrifuge at 14,000 rpm for 2 min.
 3. Remove and discard supernatant.
 4. Add 0.5 ml of Denaturing Equilibration Buffer (pH 8.0).
 5. Vortex until cell pellet is completely dissolved.
 6. Centrifuge at 14,000 rpm for 5 min to pellet any insoluble debris.
 7. Set aside 50 μ l of the supernatant for later analysis. Transfer the remainder of the supernatant to a clean 1.5-ml tube containing 50 μ l of prewashed BD TALON™ Resin, prepared as described in Section VII.B. Steps 1–7. Start with 100 μ l of resuspended slurry.
 8. Agitate sample at room temperature for 10 min.
 9. Centrifuge at 14,000 rpm for 1 min to pellet protein/resin complexes.
 10. Carefully remove the supernatant and set aside 50 μ l for later analysis. A high protein concentration in this sample indicates a problem with protein binding.
 11. Add 1 ml of Denaturing Equilibration Buffer.
 12. Vortex for a few seconds.
 13. Centrifuge at 14,000 rpm for 1 min to pellet resin.
 14. Remove the supernatant and set aside 50 μ l (“first wash”) for later analysis. Discard the remainder of the supernatant.
 15. Repeat Steps 11–14. Set aside 50 μ l for analysis.

Appendix B. Mini-Scale Protein Purification...*continued*

16. Elute bound polyhistidine-tagged protein by adding 50 μ l of Elution Buffer to the resin/protein pellet and briefly vortexing.
17. Centrifuge briefly at 14,000 rpm.
18. Carefully remove the supernatant containing the polyhistidine-tagged protein.
19. Repeat the Steps 16–18. Alternatively, if you only intend to determine the concentration of polyhistidine-tagged protein in your sample, you can achieve a more complete elution, and thus, a more accurate protein quantification by eluting with EDTA as follows:
 - a. Add 50 μ l of 100 mM EDTA (pH 8.0) and vortex briefly.
 - b. Centrifuge briefly at 14,000 rpm.
 - c. Carefully remove the supernatant containing the 6xHis protein.

Note: EDTA removes bound metal from the resin: the protein sample will contain cobalt, and the BD TALON™ Resin cannot be reused.
20. Add 12 μ l of 5X SDS/PAGE sample buffer to each of the saved samples.

Note: The sample buffer will reduce multimers to monomers; thus, only a single band will be visible on an SDS-PAGE gel, even for naturally homologous multimeric proteins.
21. Heat samples at 95–98°C for 5 min.
22. Load samples and analyze on an SDS/PAGE gel.

Appendix C: Vector Information

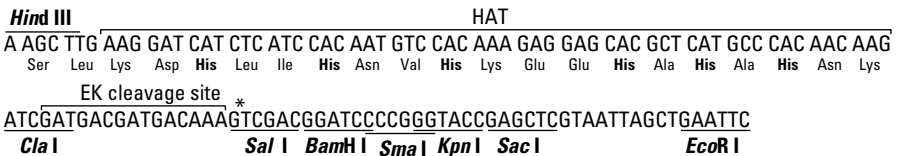
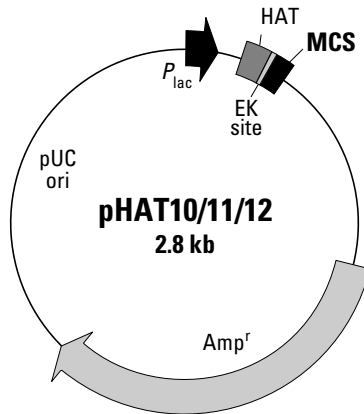


Figure 5. pHAT10/11/12 Combined Vector Map and MCS. Unique restriction sites are in bold. The sequence of pHAT10 is shown. The asterisk indicates the insertion point of additional bases in pHAT11 (G) and pHAT12 (GG) that alter the reading frame of the MCS. These vectors encode a novel polyhistidine epitope tag that enables purification of expressed proteins at neutral pH. The pHAT Vectors allow protein purification under both native and denaturing conditions. The HAT epitope is a naturally occurring, 19-amino-acid sequence from the chicken lactate dehydrogenase protein. This sequence of nonadjacent histidine residues has lower overall charge than tags with consecutive His residues, such as the 6xHis tag. As a result, HAT-protein fusions exhibit solubility that more closely resembles wild-type proteins while still possessing strong affinity for immobilized metal ions. The unique binding characteristics of the HAT sequence allow both imidazole- and pH-gradient purification of proteins under native conditions at neutral pH (7.0), as well as under denaturing conditions. The HAT sequence and an enterokinase (EK) cleavage site have been incorporated into the pUC19 backbone. The EK site allows for optional removal of the HAT sequence from the purified protein by treatment with enterokinase. Restriction sites allow excision of the HAT sequence, with or without the EK site, for cloning in other vectors.

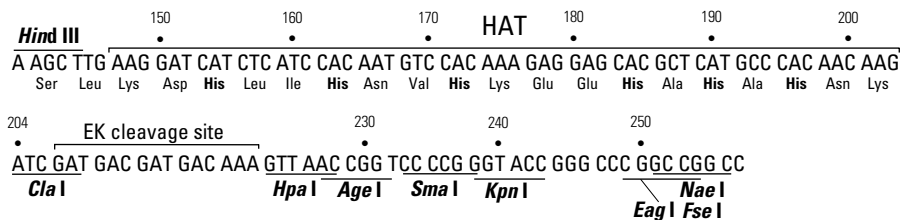
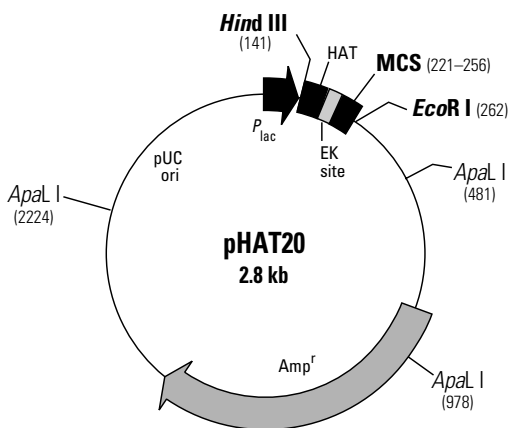
Appendix C: Vector Information *continued*

Figure 6. pHAT20 Combined Vector Map and MCS. Unique restriction sites are in bold. The sequence of pHAT20 is shown. This vector encodes a novel Histidine Affinity Tag (HAT) that enables purification of expressed proteins at neutral pH. The pHAT Vectors allow protein purification under both native and denaturing conditions. The HAT epitope is a naturally occurring, 19-amino-acid sequence from the chicken lactate dehydrogenase protein. This sequence of nonadjacent histidine residues has lower overall charge than tags with consecutive His residues, such as the 6xHis tag. As a result, HAT-protein fusions exhibit solubility that more closely resembles wild-type proteins while still possessing strong affinity for immobilized metal ions. The unique binding characteristics of the HAT sequence allow both imidazole- and pH-gradient purification of proteins under native conditions at neutral pH (7.0), as well as under denaturing conditions. The HAT sequence and an enterokinase (EK) cleavage site have been incorporated into the pUC19 backbone. The EK site allows for optional removal of the HAT sequence from the purified protein by treatment with enterokinase. Restriction sites allow excision of the HAT sequence, with or without the EK site, for cloning in other vectors.

Notes
